

# *From llama to nanobody: a streamlined workflow for the generation of functionalised VHHs*

Article

Published Version

Creative Commons: Attribution-Noncommercial 4.0

Open Access

Eyssen, L. E. -A., Ramadurai, S., Abdelkarim, S., Buckle, I., Cornish, K., Lin, H., Jones, A. K., Stephens, G. J. ORCID: <https://orcid.org/0000-0002-8966-4238> and Owens, R. J. (2024) From llama to nanobody: a streamlined workflow for the generation of functionalised VHHs. Bio-Protocol, 14 (6). e4962. ISSN 2331-8325 doi: 10.21769/BioProtoc.4962 Available at <https://centaur.reading.ac.uk/115863/>

It is advisable to refer to the publisher's version if you intend to cite from the work. See [Guidance on citing](#).

Identification Number/DOI: 10.21769/BioProtoc.4962  
<<https://doi.org/10.21769/BioProtoc.4962>>

Publisher: Bio-protocol LLC

All outputs in CentAUR are protected by Intellectual Property Rights law, including copyright law. Copyright and IPR is retained by the creators or other copyright holders. Terms and conditions for use of this material are defined in the [End User Agreement](#).

[www.reading.ac.uk/centaur](http://www.reading.ac.uk/centaur)

## **CentAUR**

Central Archive at the University of Reading

Reading's research outputs online

# From Llama to Nanobody: A Streamlined Workflow for the Generation of Functionalised VHHs

Lauren E.-A. Eyssen<sup>1,\*</sup>, Siva Ramadurai<sup>1</sup>, Sahar Abdelkarim<sup>1, Sa</sup>, Imogen Buckle<sup>1</sup>, Katy Cornish<sup>1, Sb</sup>, Hong Lin<sup>2</sup>, A.K. Jones<sup>3</sup>, Gary J. Stephens<sup>2</sup> and Raymond J. Owens<sup>1, 4,\*</sup>

<sup>1</sup>Structural Biology, Rosalind Franklin Institute, Didcot, UK

<sup>2</sup>School of Pharmacy, University of Reading, Reading, UK

<sup>3</sup>Centre for Dairy Research, University of Reading, Reading, UK

<sup>4</sup>Division of Structural Biology, Nuffield Department of Medicine, University of Oxford, Oxford, UK

<sup>Sa</sup>Department of Biochemistry, University of Oxford, Oxford, UK

<sup>Sb</sup>Department of Chemistry, University of York, York, UK

\*For correspondence: [lauren.eyssen@rffi.ac.uk](mailto:lauren.eyssen@rffi.ac.uk); [ray.owens@strubi.ox.ac.uk](mailto:ray.owens@strubi.ox.ac.uk)

## Abstract


Nanobodies are recombinant antigen-specific single domain antibodies (VHHs) derived from the heavy chain-only subset of camelid immunoglobulins. Their small molecular size, facile expression, high affinity, and stability have combined to make them unique targeting reagents with numerous applications in the biomedical sciences. From our work in producing nanobodies to over sixty different proteins, we present a standardised workflow for nanobody discovery from llama immunisation, library building, panning, and small-scale expression for prioritisation of binding clones. In addition, we introduce our suites of mammalian and bacterial vectors, which can be used to functionalise selected nanobodies for various applications such as in imaging and purification.

## Key features

- Standardise the process of building nanobody libraries and finding nanobody binders so that it can be repeated in any lab with reasonable equipment.
- Introduce two suites of vectors to functionalise nanobodies for production in either bacterial or mammalian cells.

**Keywords:** Camelid heavy chain-only antibody, Nanobody, Phage display, Nanobody labelling

## Graphical overview

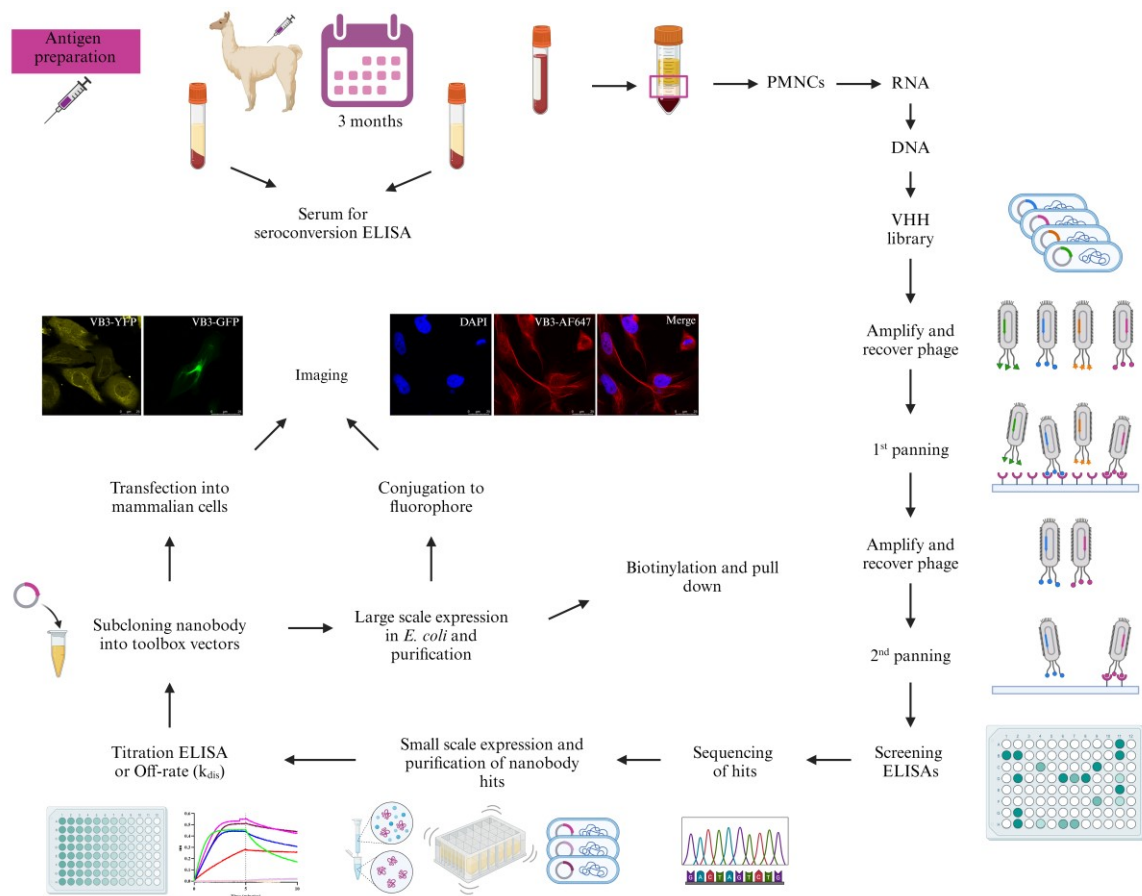

## Background

The fact that camelids produce a unique heavy-chain antibody was discovered serendipitously some thirty years ago by researchers at the Vrije Universiteit Brussel [1]. Pioneering work from the Belgian group showed that the variable heavy chain domain of these antibodies, termed VHH, could be produced as a single domain binding protein, referred to as a nanobody. Subsequently, nanobodies have been generated to a wide variety of antigens for applications in cell and structural biology, including as crystallization chaperones for high-value membrane and unstructured proteins and as probes for super-resolution microscopy [2]. Typically, nanobodies are generated by screening phage display libraries of VHH domains cloned from the peripheral blood cells of immunised llamas and alpacas with the target immunogen [3]. Since 2019, our group have generated nanobodies to 75 different antigens. The antigens varied from complexes, membranes, and soluble proteins. Some of the nanobodies that we have identified have been applied as structural chaperones [4], diagnostics [5], and anti-viral therapeutics for SARS-CoV-2 [6] and for in vitro and in vivo cell biology [7]. Although several protocols for producing nanobodies have been published previously [8], we have identified areas for streamlining the process, e.g., incorporating ligation-independent cloning to facilitate the construction of VHH domain libraries and small-scale expression screening to identify high-producing clones.

Furthermore, by using generic VHH cloning primers, we have designed a suite of expression vectors that enable the functionalisation of any cloned VHH with a variety of carboxy terminal tags, e.g., Flag, Avi-tag<sup>®</sup>, or SNAP-tag<sup>®</sup> for subsequent site-specific labelling. Overall, our aim is to make this technology readily accessible to any research group with an appropriately equipped laboratory.

## Materials and reagents

### Biological materials

1. Post-immune heparinised llama blood
2. Pre- and post-immune llama sera
3. Electrocompetent TG1 (Agilent, catalog number: 200123)
4. CM13K trypsin-sensitive helper phages (Antibody Design Laboratories, catalog number: PH050L)
5. Stellar<sup>TM</sup> competent cells (Takara Bio, catalog number: 636766)
6. *Escherichia coli* (Migula) Castellani and Chalmers, strain WK6 (ATCC, catalog number: 47078)
7. BL21(DE3)-R3-pRARE2-BirA *E. coli* cells for in vivo biotinylation, SGC (Structural Genomics Consortium)
8. HeLa cells (ATCC, catalog number: CCL-2)

### Reagents

1. Gerbu adjuvant F (Biotechnik GmbH, catalog number: 3030). Store at 4 °C
2. NeutrAvidin<sup>TM</sup> biotin-binding protein (Invitrogen, catalog number: A2666). Store at -20 °C
3. 10× PBS (Fisher BioReagents, catalog number: BP399-20), for general use. Store at room temperature (RT)
4. Skim milk powder (Oxoid, catalog number: LP0031). Store at RT
5. Tween-20 (Sigma-Aldrich, catalog number: P1379). Store at RT
6. BSA (Sigma-Aldrich, catalog number: A2153). Store at 4 °C
7. Goat anti-llama IgG (H+L) HRP (Invitrogen, catalog number: A16060). Store at -20 °C in 5 µL aliquots
8. KPL ABTS peroxidase solution A (SeraCare, catalog number: 5120-0034). Store at 4 °C
9. KPL peroxidase substrate solution B (SeraCare, catalog number: 5120-0037). Store at 4 °C
10. Chemgene<sup>TM</sup> (Chemgene&Trade, catalog number: XTM308). Store at RT  
**Caution:** This chemical is corrosive.
11. Ethanol (EtOH) (Fisher BioReagents, catalog number: E/0650DF/17). Store at RT in a flammable liquid storage cupboard  
**Caution:** This chemical is flammable.
12. Virkon<sup>®</sup> (Rely+on<sup>TM</sup>, catalog number: 12358667). Store at RT  
**Caution:** This chemical is corrosive.
13. PBS, pH 7.4 for polymorphonuclear cell (PMNC) isolation (Gibco, catalog number: 10010015). Store at 4 °C
14. Histopaque<sup>®</sup>-1077 (Sigma-Aldrich, catalog number: 10771). Store at 4 °C
15. RNaseZap<sup>TM</sup> RNase decontamination solution (Invitrogen, catalog number: AM9782). Store at RT
16. Trypan blue solution, 0.4% (Sigma-Aldrich, catalog number: T8154-20ML). Store at RT  
**Caution:** This chemical is potentially carcinogenic.
17. TRIzol<sup>TM</sup> reagent (Invitrogen, catalog number: 15596026). Store at 4 °C  
**Caution:** This chemical is corrosive and toxic.
18. Chloroform (Fluorchem, catalog number: D007F). Store at RT in a flammable liquid storage cupboard  
**Caution:** This chemical is flammable.
19. 2-propanol (Sigma-Aldrich, catalog number: I9516). Store at RT in a flammable liquid storage cupboard  
**Caution:** This chemical is flammable.
20. Nuclease-free water (not DEPC treated) (Invitrogen, catalog number: AM9932). Store at RT
21. SuperScript<sup>TM</sup> IV one-step RT-PCR system (Invitrogen, catalog number: 12594025). Store at -20 °C
22. CALL\_001 primer 5' GTCTGGCTGCTCTTCTACAAGG 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
23. CALL\_002 primer 5' GGTACGTGCTGTTGAAGTGTCC 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
24. 6× DNA gel loading buffer (New England Biolabs, catalog number: B7024S). Store at 4 °C
25. 10× TBE buffer (Thermo Scientific, catalog number: B52). Store at RT
26. Agarose (Fisher BioReagents, catalog number: BP1356-500). Store at RT

27. 1,000× SYBR<sup>TM</sup> safe DNA gel stain (Invitrogen, catalog number: S33102). Store at RT
28. GeneRuler 1 kb DNA ladder (Thermo Scientific, catalog number: SM0311). Store at 4 °C
29. Nucleospin Gel and PCR clean-up kit (Macherey-Nagel, catalog number: 740609.50). Store at RT
30. Purelink PCR purification kit (Invitrogen, catalog number: K310002). Store at RT
31. HyperLadder<sup>TM</sup> 1 kb DNA ladder (Meridian Bioscience<sup>®</sup>, Bionline, catalog number: BIO-33053). Store at 4 °C
32. Phusion flash high fidelity PCR master mix (Invitrogen, catalog number: F548L). Store at -20 °C
33. VHHFor2 primer 5' TACTCGCGGCCAGCCGGCCATGGCCAGGTGCAGCTGCAGGAGTCT GGRGGA 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
34. VHHRev primer 5' GTGATGGTGTGGCCTCCTGAGGAGACGGTGACCTGG 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
35. pADL23c vector (Antibody Design Laboratories, catalog number: PD0111). Store at -20 °C
36. SfiI (New England Biolabs, catalog number: R0123S). Store at -20 °C
37. 10× rCutSmart<sup>TM</sup> buffer (New England Biolabs, catalog number: B6004S). Store at 4 °C
38. ClonExpress II one step cloning kit, which contains 5× CEII buffer and Exnase II (Vazymbiotech, catalog number: C112-02). Store at -20 °C
39. Recovery medium (Sigma-Aldrich, catalog number: S1797). Store at 4 °C
40. LB medium mix (Formedium, catalog number: LBL0103). Store at RT
41. Bacto tryptone (Melford, catalog number: T60065-2000.0). Store at RT
42. Yeast extract (Melford, catalog number: Y20020-1000.0). Store at RT
43. Bacto agar (Formedium, catalog number: AGR10). Store at RT
44. Ampicillin (Formedium, catalog number: AMP100). Store at 4 °C  
**Caution:** Ampicillin is a sensitiser.
45. 25 mM dNTP (Thermo Scientific, catalog number: R1122). Store at -20 °C in 30 µL aliquots
46. *Taq* polymerase with standard *Taq* reaction buffer (New England Biolabs, catalog number: M0273X). Store at -20 °C
47. PhD\_seq\_Fwd primer 5' GCTTCCGGCTCGTATGTTG 3'. Store at -20 °C
48. PhD\_seq\_Rev primer 5' GTCGTCTTTCCAGACGTTAG 3'. Store at -20 °C
49. Glycerol (Sigma-Aldrich, catalog number: G5516-1L). Store at RT
50. PEG 6000 (Sigma-Aldrich, catalog number: 81260). Store at RT
51. NaCl (Sigma-Aldrich, catalog number: S9888). Store at RT
52. Kanamycin (Sigma-Aldrich, catalog number: K1377). Store at 4 °C  
**Caution:** Kanamycin is a sensitiser and may damage fertility.
53. StartingBlock<sup>TM</sup> (PBS) blocking buffer (Invitrogen, catalog number: 37538). Store at 4 °C
54. Dynabeads<sup>TM</sup> M-280 streptavidin (Invitrogen, catalog number: 11205D). Store at 4 °C
55. Tris base (Melford, catalog number: T60040-100.0). Store at RT
56. CaCl<sub>2</sub>·2H<sub>2</sub>O (Sigma-Aldrich, catalog number: C3306). Store at RT
57. Trypsin (Sigma-Aldrich, catalog number: T1426). Store at -20 °C
58. Anti-M13 HRP (Sino Biological, catalog number: 11973-MM05T-H). Store at -20 °C in 5 µL aliquots
59. Highprep PCR magnetic beads (Magbio, catalog number: AC60050). Store at 4 °C
60. Glucose (Sigma-Aldrich, catalog number: G8270). Store at RT
61. MgCl<sub>2</sub>·6H<sub>2</sub>O (Sigma-Aldrich, catalog number: M2670). Store at RT
62. IPTG (NeoBiotech, catalog number: NB-45-00030). Store at -20 °C
63. Polymyxin B sulfate (Gibco, catalog number: 21850029). Store at RT  
**Caution:** Polymyxin B sulfate is toxic.
64. 2× Laemmli buffer (Sigma-Aldrich, catalog number: S3401). Store at 4 °C
65. 20× NuPAGE<sup>TM</sup> MES SDS running buffer (Invitrogen, catalog number: NP000202). Store at RT
66. Mark12<sup>TM</sup> unstained standard (Invitrogen, catalog number: LC5677). Store at 4 °C
67. InstantBlue<sup>®</sup> Coomassie protein stain (Abcam, catalog number: ab119211). Store at 4 °C
68. Ni-NTA spin columns (QIAGEN, catalog number: 31014). Store at 4 °C
69. Imidazole (Sigma-Aldrich, catalog number: I202-500g). Store at RT  
**Caution:** Imidazole is toxic, corrosive, and an irritant and may damage fertility.
70. Zeba<sup>TM</sup> spin desalting columns, 7K MWCO, 0.5 mL (Thermo Scientific, catalog number: 89882). Store at 4 °C

71. Rabbit anti-camelid VHH HRP (GenScript, catalog number: A01861-200). Store at -20 °C in 5 µL aliquots
72. pOPINVHH\_his vector (Addgene, catalog number: 210405). Store at -20 °C
73. pOPINVHH\_cys\_his vector (Addgene, catalog number: 210403). Store at -20 °C
74. pOPINVHH\_BAP\_his vector (Addgene, catalog number: 210402). Store at -20 °C
75. pOPINVHH\_sort\_his vector (Addgene, catalog number: 210406). Store at -20 °C
76. pOPINVHH\_flag\_his vector (Addgene, catalog number: 210404). Store at -20 °C
77. pOPINVHH\_myc\_his vector (Addgene, catalog number: 210526). Store at -20 °C
78. pOPINVHH\_snap\_his vector (Addgene, catalog number: 210407). Store at -20 °C
79. pOPINE-3C-eGFP vector (Addgene, catalog number: 41125). Store at -20 °C
80. pOPINE-3C-eYFP vector (Addgene, catalog number: 214028). Store at -20 °C
81. pOPINE-3C-mCherry vector (Addgene, catalog number: 214060). Store at -20 °C
82. Common Fwd primer 5' GCGGCCAGCCGCCATGCCCCAGGTGCAGCTGGTGGAG 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
83. His FLAG Rev primer 5' GTGATGGTGGCCTGAGGAGACGGTGACCTGGGTC 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
84. Cys Rev primer 5' ATGGTGACAGCCTGAGGAGACGGTGACCTGGGTC 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
85. BAP Rev primer 5' ATCATTCAAGCCTGAGGAGACGGTGACCTGGGTC 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C.
86. Sort Rev primer 5' CGGCAGGCCGCCTGAGGAGACGGTGACCTGGGTC 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
87. Myc Rev primer 5' GTGATGGTGGTGGCCTGAGGAGACGGTGACCTGGGTC 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
88. SNAP Rev primer 5' GTCCTTGTCGCCTGAGGAGACGGTGACCTGGGTC 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
89. pOPINE common Fwd primer 5' AGGAGATATACCATGCAGGTGCAGCTGGTGGAG 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
90. pOPINE common Rev primer 5' CAGAACTTCCAGTTTAGGGGAGACGGTGACCTGGGTC 3'. Ordered primer from IDT. Store at RT prior to reconstitution in water, then store at -20 °C
91. Bsu36I (New England Biolabs, catalog number: R0524S). Store at -20 °C
92. NcoI (New England Biolabs, catalog number: R01093S). Store at -20 °C
93. X-gal ready to use (Thermo Scientific, catalog number: R0941). Store at 4 °C
94. QIAprep spin miniprep kit (QIAGEN, catalog number: 27104). Store at RT
95. DNA sequence of anti-GFP nanobody [9] after codon optimisation for expression in *E. coli*. CAG GTC CAA TTA GTG GAG TCC GGT GGG GCA CTT GTC CAG CCT GGA GGT TCA CTT CGC TTG TCT TGC GCA GCG TCT GGA TTC CCG GTG AAC CGC TAT AGT ATG CGT TGG TAC CGT CAA GCT CCG GGG AAA GAA CGT GAA TGG GTA GCA GGG ATG TCT TCC GCC GGT GAC CGC TCT TCA TAC GAG GAC TCG GTC AAG GGG CGC TTC ACA ATC TCT CGT GAT GAT GCC CGT AAC ACC GTT TAC TTG CAA ATG AAC AGC CTG AAA CCG GAA GAC ACT GCG GTG TAT TAC TGC AAT GTT AAT GTA GGG TTT GAA TAC TGG GGT CAA GGT ACA CAA GTT ACA GTT TCG TCA. Ordered gBlock from IDT. Store at RT prior to reconstitution in TE buffer, incubate at 50 °C for 20 min, then store at -20 °C
96. Chloramphenicol (Sigma-Aldrich, catalog number: C1919). Store at 4 °C  
**Caution:** Chloramphenicol is toxic, can cause eye damage, and is a carcinogen.
97. Spectinomycin dihydrochloride pentahydrate (Sigma-Aldrich, catalog number: S4014). Store at 4 °C  
**Caution:** Spectinomycin is an irritant.
98. Carbenicillin disodium salt (Sigma-Aldrich, catalog number: C3416). Store at 4 °C  
**Caution:** Carbenicillin is a sensitizer.
99. Biotin (Sigma-Aldrich, catalog number: B4639). Store at 4 °C
100. Streptavidin, Alexa Fluor™ 488 conjugate (Invitrogen, catalog number: S11223). Store at -20 °C
101. BenchMark™ fluorescent protein standard (Invitrogen, catalog number: LC5928). Store at -20 °C
102. DNA sequence of anti-vimentin (VB3) nanobody [10] after codon optimisation for expression in *E. coli*: CAG GTC CAA CTT GTA GAG TCA GGA GGT GGA AGC GTG CAA GCT GGG GAC TCT CTG CGC CTG



TCT TGT GCT TCG AGC GGA AAT ACC TTC TCG ATC AAA GTC ATG GGA TGG TAC CGC CAG GCA CCT GGA AAG CAA CGT GAA TTA GTC GCG GTT TCA ACC AAT AGC GGG GCC TCT GTT AAT TAT GCC AAC TCT GTG AAG GGA CGC TTT ACC ATT TCT ATT GAT TCA GTA AAA AAA ACA ACC TAC TTA CAG ATG AAT TCC TTG AAG CCA GAA GAT ACA GCC GTC TAC TTT TGC AAT GCA TAT GAT GGG CGT TAT GAG GAC TAT TAC GGT CAG GGG ACC CAA GTG ACA GTA TCA TCA. Ordered gBlock from IDT. Store at RT prior to reconstitution in TE buffer, incubate at 50 °C for 20 min, then store at -20 °C

103. TCEP, hydrochloride (Merck, Millipore<sup>®</sup>, catalog number: 580567-5GM). Store at RT
104. Alexa Fluor<sup>™</sup> 647 C<sub>2</sub>-maleimide (Invitrogen, catalog number: A20347). Store at -20 °C
105. Zeba<sup>™</sup> dye and biotin removal columns (Thermo Scientific, catalog number: A44296). Store at 4 °C
106. PageRuler<sup>™</sup> prestained protein ladder (Thermo Scientific, catalog number: 26616). Store at 4 °C
107. DMEM, high glucose, HEPES, no phenol red (Gibco, catalog number: 21063-029). Store at 4 °C
108. Fetal bovine serum (FBS) (Biowest, catalog number: S00NB1001Y). Store at -20 °C
109. 100× GlutaMAX<sup>™</sup> (Gibco, catalog number: 35050-038). Store at -20 °C
110. 100× penicillin-streptomycin (Gibco, catalog number: 15140-122). Store at -20 °C
111. 4% paraformaldehyde in PBS (Thermo Scientific chemicals, catalog number: J61899.AK). Store at 4 °C
- Caution:** Paraformaldehyde is an eye irritant, skin sensitizer, and carcinogen.
112. Triton<sup>™</sup> X-100 (Sigma-Aldrich, catalog number: T8787-250mL). Store at RT
113. Fluoroshield with DAPI (Abcam, catalog number: ab104139). Store at RT
114. Nail varnish. Store at RT
115. Type F immersion liquid (Leica, catalog number: 11513859). Store at RT
116. Eugene transfection reagent (Promega, catalog number: E2311). Store at 4 °C

## Solutions

1. 0.5 mg/mL neutrAvidin biotin-binding protein (see Recipes)
2. 10 µg/mL neutrAvidin for coating ELISA plate (see Recipes)
3. 1× PBS (see Recipes)
4. Blocking solution (see Recipes)
5. Washing buffer (PBST) (see Recipes)
6. 0.1% (w/v) BSA-PBS (see Recipes)
7. 5% (v/v) Chemgene (see Recipes)
8. 70% (v/v) ethanol (EtOH) (see Recipes)
9. 2% (w/v) Virkon (see Recipes)
10. 75% (v/v) EtOH (see Recipes)
11. 1× TBE buffer (see Recipes)
12. 0.7% (w/v) agarose gel containing 1× SYBR<sup>™</sup> Safe DNA gel stain (see Recipes)
13. 1% (w/v) agarose gel containing 1× SYBR<sup>™</sup> Safe DNA gel stain (see Recipes)
14. 2× YT medium (see Recipes)
15. LB medium (see Recipes)
16. 100 mg/mL ampicillin (see Recipes)
17. 1% (w/v) agar LB plates containing 100 µg/mL ampicillin (see Recipes)
18. 2% (w/v) agar LB plates containing 100 µg/mL ampicillin (see Recipes)
19. 50% glycerol (see Recipes)
20. 2× YT medium containing 25% glycerol (see Recipes)
21. 1% (w/v) agar LB plates without antibiotic (see Recipes)
22. 2× YT containing 100 µg/mL ampicillin (see Recipes)
23. 50 mg/mL kanamycin (see Recipes)
24. 2× YT containing 25 µg/mL kanamycin (see Recipes)
25. PEG/NaCl precipitation solution (see Recipes)
26. 2× YT containing 100 µg/mL ampicillin and 25 µg/mL kanamycin (see Recipes)
27. TBSC (see Recipes)
28. 1 mg/mL trypsin (see Recipes)



29. 250 µg/mL trypsin (see Recipes)
30. Terrific broth (see Recipes)
31. 20% (w/v) glucose (see Recipes)
32. 1 M MgCl<sub>2</sub>·6H<sub>2</sub>O (see Recipes)
33. Terrific broth containing 100 µg/mL ampicillin, 0.1% glucose, 2 mM MgCl<sub>2</sub>·6H<sub>2</sub>O (see Recipes)
34. 1 M IPTG (see Recipes)
35. 1 mg/mL polymyxin B sulfate in PBS (see Recipes)
36. 1× NuPAGE™ MES SDS running buffer (see Recipes)
37. Equilibration buffer (see Recipes)
38. Elution buffer (see Recipes)
39. Octet dilution buffer (see Recipes)
40. 1% (w/v) agar LB plate containing 100 µg/mL ampicillin, 2 mM IPTG, 40 µg/mL X-gal (see Recipes)
41. 34 µg/mL chloramphenicol (see Recipes)
42. 50 µg/mL spectinomycin (see Recipes)
43. 50 µg/mL carbenicillin (see Recipes)
44. 1% (w/v) agar LB plate containing 34 µg/mL chloramphenicol, 50 µg/mL spectinomycin, 50 µg/mL carbenicillin (see Recipes)
45. 0.2 M Biotin (see Recipes)
46. Terrific broth containing 50 µg/mL spectinomycin and 50 µg/mL carbenicillin (see Recipes)
47. 1 mg/mL Streptavidin Alexa Fluor™ 488 conjugate (see Recipes)
48. 50 µg/mL Streptavidin Alexa Fluor™ 488 conjugate (see Recipes)
49. 500 mM TCEP, pH 8.0 (see Recipes)
50. 1 mM TCEP, pH 8.0 (see Recipes)
51. 10 mg/mL Alexa Fluor™ 647 C<sub>2</sub>-maleimide (see Recipes)
52. DMEM medium without phenol red, 10% (v/v) FBS, 1× GlutaMAX™, 1× Penicillin-streptomycin (see Recipes)
53. Confocal blocking buffer (see Recipes)
54. Confocal dilution buffer (see Recipes)

## Recipes

### 1. 0.5 mg/mL neutrAvidin biotin-binding protein

Prepare 200 µL aliquots and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
neutrAvidin™	0.5 mg	5 mg
ddH <sub>2</sub> O	n/a	10 mL
Total	n/a	10 mL

### 2. 10 µg/mL neutrAvidin for coating ELISA plate

Prepare just before use.

Reagent	Final concentration	Quantity or Volume
neutrAvidin™ (Recipe 1)	10 µg/mL	200 µL
PBS (1×)	1×	9.8 mL
Total	n/a	10 mL

### 3. 1× PBS

Prepare both sterile (by autoclaving) and non-sterile PBS. Store at RT. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
PBS (10×)	1×	100 mL
ddH <sub>2</sub> O	n/a	900 mL
Total	n/a	1,000 mL

**4. Blocking solution**

Prepare on day of use.

Reagent	Final concentration	Quantity or Volume
Milk	2% (w/v)	2 g
PBS (Recipe 3) (or buffer that is compatible with the protein)	n/a	100 mL
Total	n/a	100 mL

**5. Washing buffer (PBST)**

Store at RT. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
PBS (Recipe 3)	1×	1 000 mL
Tween-20	0.05% (v/v)	500 µL
Total	n/a	1,000 mL

**6. 0.1% (w/v) BSA-PBS**

Prepare on day of use.

Reagent	Final concentration	Quantity or Volume
BSA	0.1 % (w/v)	50 mg
PBS (Recipe 3)	1×	50 mL
Total	n/a	50 mL

**7. 5% (v/v) Chemgene**

Store at RT.

Reagent	Final concentration	Quantity or Volume
Chemgene (100%)	5% (v/v)	50 mL
H <sub>2</sub> O	n/a	950 mL
Total	n/a	1,000 mL

**8. 70% (v/v) EtOH**

Store at RT.

Reagent	Final concentration	Quantity or Volume
Ethanol (absolute)	70% (v/v)	70 mL
H <sub>2</sub> O	n/a	30 mL
Total	n/a	100 mL

**9. 2% (w/v) Virkon**

Prepare weekly and store at RT.

Reagent	Final concentration	Quantity or Volume
Virkon powder	2% (w/v)	100 g
H <sub>2</sub> O	n/a	5,000 mL
Total	n/a	5,000 mL

**10. 75% (v/v) EtOH**

Store at RT.

Reagent	Final concentration	Quantity or Volume
Ethanol (absolute)	75% (v/v)	75 mL
H <sub>2</sub> O	n/a	25 mL
Total	n/a	100 mL

#### 11. 1× TBE buffer

Store at RT. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
TBE (10×)	1×	100 mL
H <sub>2</sub> O	n/a	900 mL
Total	n/a	1,000 mL

#### 12. 0.7% (w/v) agarose gel containing 1× SYBR<sup>TM</sup> Safe DNA gel stain

Heat in a microwave until all the agarose has dissolved and allow to cool to 55 °C before the addition of the gel stain. Pour into the DNA casting apparatus and allow to solidify. Prepare on day of use.

Reagent	Final concentration	Quantity or Volume
Agarose	0.7 % (w/v)	0.7 g
TBE (Recipe 11)	1×	100 mL
Total	n/a	100 mL
SYBR <sup>TM</sup> Safe DNA gel stain (1,000×)	1×	10 µL

#### 13. 1% (w/v) agarose gel containing 1× SYBR<sup>TM</sup> Safe DNA gel stain

Heat in a microwave until all the agarose has dissolved and allow to cool to 55 °C before the addition of the gel stain. Pour into the DNA casting apparatus and allow to solidify. Prepare on day of use.

Reagent	Final concentration	Quantity or Volume
Agarose	1% (w/v)	1 g
TBE (Recipe 11)	1×	100 mL
Total	n/a	100 mL
SYBR <sup>TM</sup> Safe DNA gel stain (1,000×)	1×	10 µL

#### 14. 2× YT medium

Autoclave and store at RT. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
Bacto tryptone	16 g/L	4 g
NaCl	5 g/L	1.25 g
Yeast extract	10 g/L	2.5 g
ddH <sub>2</sub> O	n/a	250 mL
Total	n/a	250 mL

#### 15. LB medium

Autoclave and store at RT. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
LB medium mix	25 g/L	6.25 g
ddH <sub>2</sub> O	n/a	250 mL
Total	n/a	250 mL

#### 16. 100 mg/mL ampicillin

Filter sterilise using a 0.22 µm filter. Prepare 500 µL aliquots and store at -20 °C. Stable for at least one year.

**Caution:** Ampicillin is a sensitizer.

Reagent	Final concentration	Quantity or Volume
Ampicillin	100 mg/mL	2 g
ddH <sub>2</sub> O	n/a	20 mL
Total	n/a	20 mL

#### 17. 1% (w/v) agar LB plates containing 100 µg/mL ampicillin

Autoclave LB medium and bacto agar mix and allow to cool to 55 °C before the addition of antibiotic. Invert

to mix and pour into 8.5 cm Petri dishes in a safety cabinet. Allow to cool at RT. Prepare on day of use.

Reagent	Final concentration	Quantity or Volume
LB medium mix	25 g/L	6.25 g
Bacto agar	1% (w/v)	2.5 g
ddH <sub>2</sub> O	n/a	250 mL
Total	n/a	250 mL
Ampicillin (Recipe 16)	100 µg/mL	250 µL

#### 18. 2% (w/v) agar LB plates with 100 µg/mL ampicillin

Autoclave LB medium and bacto agar mix and allow to cool to 55 °C before the addition of antibiotic. Invert to mix and pour into 8.5 cm Petri dishes in a safety cabinet. Allow to cool at RT. Prepare on day of use.

Reagent	Final concentration	Quantity or Volume
LB medium mix	25 g/L	6.25 g
Bacto agar	2% (w/v)	5 g
ddH <sub>2</sub> O	n/a	250 mL
Total	n/a	250 mL
Ampicillin (Recipe 16)	100 µg/mL	250 µL

#### 19. 50% (v/v) glycerol

Autoclave and store at RT. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
Glycerol	50% (v/v)	100 mL
ddH <sub>2</sub> O	n/a	100 mL
Total	n/a	200 mL

#### 20. 2× YT medium containing 25% (v/v) glycerol

Prepare just before use in a safety cabinet.

Reagent	Final concentration	Quantity or Volume
2× YT medium	n/a	20 mL
Glycerol (Recipe 19)	25% (w/v)	20 mL
Total	n/a	40 mL

#### 21. 1% (w/v) agar LB plates without antibiotic

Autoclave and allow to cool to 55 °C. Invert to mix and pour into 8.5 cm Petri dishes in a safety cabinet. Allow to cool and set at RT. Store at 4 °C for two weeks.

Reagent	Final concentration	Quantity or Volume
LB medium mix	25 g/L	6.25 g
Bacto agar	1% (w/v)	2.5 g
ddH <sub>2</sub> O	n/a	250 mL
Total	n/a	250 mL

#### 22. 2× YT containing 100 µg/mL ampicillin

Prepare just before use in a safety cabinet.

Reagent	Final concentration	Quantity or Volume
2× YT	n/a	100 mL
Ampicillin (Recipe 16)	100 µg/mL	100 µL
Total	n/a	100 mL

#### 23. 50 mg/mL kanamycin

Filter sterilise using a 0.22 µm filter. Prepare 100 µL aliquots and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
Kanamycin	50 mg/mL	0.5 g
ddH <sub>2</sub> O	n/a	10 mL
Total	n/a	10 mL

**24. 2× YT containing 25 µg/mL kanamycin**

Prepare just before use in a safety cabinet.

Reagent	Final concentration	Quantity or Volume
2× YT	n/a	100 mL
Kanamycin (Recipe 23)	25 µg/mL	50 µL
Total	n/a	100 mL

**25. PEG/NaCl precipitation solution**

Autoclave and store at RT. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
PEG 6000	20% (w/v)	200 g
NaCl	2.5 M	146.1 g
ddH <sub>2</sub> O	n/a	Make up to 1,000 mL
Total	n/a	1,000 mL

**26. 2× YT containing 100 µg/mL ampicillin and 25 µg/mL kanamycin**

Prepare just before use in a safety cabinet.

Reagent	Final concentration	Quantity or Volume
2× YT	n/a	100 mL
Ampicillin (Recipe 16)	100 µg/mL	100 µL
Kanamycin (Recipe 23)	25 µg/mL	50 µL
Total	n/a	100 mL

**27. TBSC**

Autoclave and store at RT. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
Tris base	10 mM	0.788 g
NaCl	137 mM	4 g
CaCl <sub>2</sub> ·2H <sub>2</sub> O	1 mM	73.5 mg
ddH <sub>2</sub> O	n/a	500 mL
Total	n/a	500 mL

**28. 1 mg/mL trypsin**

Prepare 130 µL aliquots and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
Trypsin	1 mg/mL	50 mg
ddH <sub>2</sub> O	n/a	50 mL
Total	n/a	50 mL

**29. 250 µg/mL trypsin**

Prepare just before use.

Reagent	Final concentration	Quantity or Volume
Trypsin (Recipe 28)	250 µg/mL	125 µL
TBSC (Recipe 27)	n/a	375 µL
Total	n/a	500 µL

### 30. Terrific broth

Autoclave and store at RT. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
Bacto tryptone	12 g/L	3 g
Yeast extract	24 g/L	6 g
Glycerol	0.4% (v/v)	1 mL
KH <sub>2</sub> PO <sub>4</sub>	0.17 M	0.5775 g
K <sub>2</sub> HPO <sub>4</sub>	0.72 M	3.135 g
ddH <sub>2</sub> O	n/a	250 mL
Total	n/a	250 mL

### 31. 20% (w/v) glucose

Filter sterilise using a 0.22 µm filter. Prepare 500 µL aliquots and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
Glucose	20% (w/v)	20 g
ddH <sub>2</sub> O	n/a	Make up to 100 mL
Total	n/a	100 mL

### 32. 1 M MgCl<sub>2</sub>·6H<sub>2</sub>O

Filter sterilise using a 0.22 µm filter. Prepare 500 µL aliquots and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
MgCl <sub>2</sub> ·6H <sub>2</sub> O	1 M	4.066 g
ddH <sub>2</sub> O	n/a	20 mL
Total	n/a	20 mL

### 33. Terrific broth containing 100 µg/mL ampicillin, 0.1% glucose, 2 mM MgCl<sub>2</sub>·6H<sub>2</sub>O

Prepare just before use in a safety cabinet.

Reagent	Final concentration	Quantity or Volume
Terrific broth (Recipe 30)	n/a	99.2 mL
Ampicillin (Recipe 16)	100 µg/mL	100 µL
Glucose (Recipe 31)	0.1% (v/v)	500 µL
MgCl <sub>2</sub> ·6H <sub>2</sub> O (Recipe 32)	2 mM	200 µL
Total	n/a	100 mL

### 34. 1 M IPTG

Filter sterilise using a 0.22 µm filter. Prepare 500 µL aliquots and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
IPTG	1 M	4.766 g
ddH <sub>2</sub> O	n/a	20 mL
Total	n/a	20 mL

### 35. 1 mg/mL polymyxin-B sulfate in PBS

Prepare 1 mL aliquots and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
Polymyxin-B sulfate	1 mg/mL	2 mg
PBS (Recipe 3)	1×	2 mL
Total	n/a	2 mL

### 36. 1× NuPAGE™ MES SDS running buffer

Store at RT. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
MES (20×)	1×	100 mL
H <sub>2</sub> O	n/a	1,900 mL
Total	n/a	2,000 mL

### 37. Equilibration buffer

Store at RT. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
Imidazole	30 mM	0.51 g
PBS (Recipe 3)	1×	250 mL
Adjust pH to 7.4		
Total	n/a	250 mL

### 38. Elution buffer

Store at RT. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
Imidazole	300 mM	2.04 g
PBS (Recipe 3)	1×	100 mL
Adjust pH to 7.4		
Total	n/a	100 mL

### 39. Octet dilution buffer

Store at RT. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
BSA	0.1% (w/v)	50 mg
PBS (Recipe 3) (or buffer that is compatible with the protein)	1×	50 mL
Total	n/a	50 mL

### 40. 1% (w/v) agar LB plate containing 100 µg/mL ampicillin, 2 mM IPTG, 40 µg/mL X-gal

Autoclave and allow to cool to 55 °C before the addition of additives. Invert to mix and pour into 8.5 cm Petri dishes in a safety cabinet. Allow to cool and set at RT. Store at 4 °C for two weeks.

Reagent	Final concentration	Quantity or Volume
LB medium mix	25 g/L	6.25 g
Bacto agar	1% (w/v)	2.5 g
ddH <sub>2</sub> O	n/a	Make up to 250 mL
Total	n/a	250 mL
Ampicillin (Recipe 16)	100 µg/mL	250 µL
IPTG (Recipe 34)	2 mM	500 µL
X-gal (20 mg/mL)	40 µg/mL	500 µL

### 41. 34 mg/mL chloramphenicol

Filter sterilise using a 0.22 µm filter. Prepare 100 µL aliquots and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
Chloramphenicol	34 mg/mL	34 mg
EtOH	100%	1 mL
Total	n/a	1 mL

### 42. 50 mg/mL Spectinomycin

Filter sterilise using a 0.22 µm filter. Prepare 100 µL aliquots and store at -20 °C. Stable for at least one year.



Reagent	Final concentration	Quantity or Volume
Spectinomycin dihydrochloride pentahydrate	50 mg/mL	50 mg
ddH <sub>2</sub> O	n/a	1 mL
Total	n/a	1 mL

**43. 50 mg/mL Carbenicillin**

Filter sterilise using a 0.22 µm filter. Prepare 100 µL aliquots and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
Carbenicillin disodium salt	50 mg/mL	50 mg
ddH <sub>2</sub> O	n/a	1 mL
Total	n/a	1 mL

**44. 1% (w/v) agar LB plate containing 34 µg/mL chloramphenicol, 50 µg/mL spectinomycin, 50 µg/mL carbenicillin**

Autoclave LB medium and bacto agar mix and allow to cool to 55 °C before the addition of additives. Invert to mix and pour into 8.5 cm Petri dishes in a safety cabinet. Allow to cool and set at RT. Store at 4 °C for two weeks.

Reagent	Final concentration	Quantity or Volume
LB medium mix	25 g/L	2.5 g
Bacto agar	1% (w/v)	1 g
ddH <sub>2</sub> O	n/a	Make up to 100 mL
Total	n/a	100 mL
Chloramphenicol (Recipe 41)	34 µg/mL	100 µL
Spectinomycin (Recipe 42)	50 µg/mL	100 µL
Carbenicillin (Recipe 43)	50 µg/mL	100 µL

**45. 0.2 M biotin**

Filter sterilise using a 0.22 µm filter. Prepare 100 µL aliquots and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
Biotin	0.2 M	0.488 g
NaOH (1 M)	n/a	Add to dissolve biotin
ddH <sub>2</sub> O	n/a	Make up to 10 mL
Total	n/a	10 mL

**46. Terrific broth containing 50 µg/mL spectinomycin and 50 µg/mL carbenicillin**

Prepare just before use in a safety cabinet.

Reagent	Final concentration	Quantity or Volume
Terrific broth (Recipe 30)	n/a	998 mL
Spectinomycin (Recipe 42)	50 µg/mL	1 mL
Carbenicillin (Recipe 43)	50 µg/mL	1 mL
Total	n/a	1 000 mL

**47. 1 mg/mL streptavidin Alexa Fluor™ 488 conjugate**

Protect from light, prepare 10 µL aliquots, and store at -20 °C. Stable for six months.

Reagent	Final concentration	Quantity or Volume
Streptavidin Alexa Fluor™ 488 conjugate	1 mg/mL	1 mg
PBS (Recipe 3)	1×	1 mL
Total	n/a	1 mL

**48. 50 µg/mL streptavidin Alexa Fluor™ 488 conjugate**

Protect from light and prepare just before use.

Reagent	Final concentration	Quantity or Volume
Streptavidin Alexa Fluor™ 488 conjugate (1 mg/mL)	50 µg/mL	1 µL
PBS (Recipe 3)	1×	19 µL
Total	n/a	20 µL

**49. 500 mM TCEP, pH 8.0**

Store at 4 °C. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
TCEP	500 mM	1.433 g
ddH <sub>2</sub> O	n/a	10 mL
Adjust to pH 8.0		
Total	n/a	10 mL

**50. 1 mM TCEP, pH 8.0**

Prepare just before use.

Reagent	Final concentration	Quantity or Volume
TCEP (Recipe 49)	1 mM	2 µL
ddH <sub>2</sub> O	n/a	998 µL
Total	n/a	1 mL

**51. 10 mg/mL Alexa Fluor™ 647 C<sub>2</sub>-Maleimide**

Protect from light and store at -20 °C. Stable for at least one year.

Reagent	Final concentration	Quantity or Volume
Alexa Fluor™ 647 C <sub>2</sub> -Maleimide	10 mg/mL	1 mg
DMSO	1×	500 µL
Total	n/a	500 µL

**52. DMEM medium without phenol red, 10% (v/v) FBS, 1× GlutaMAX™, 1× penicillin-streptomycin**

Store at 4 °C. Stable for at least one month.

Reagent	Final concentration	Quantity or Volume
DMEM medium	n/a	88 mL
FBS	10% (v/v)	10 mL
GlutaMAX™ (100×)	1×	1 mL
Penicillin-streptomycin (100×)	1×	1 mL
Total	n/a	100 mL

**53. Confocal blocking buffer**

Store at 4 °C. Use within one month.

Reagent	Final concentration	Quantity or Volume
PBS (Recipe 3)	1×	10 mL
BSA	5% (w/v)	0.5 g
Triton X-100	0.25% (v/v)	25 µL
Total	n/a	10 mL

**54. Confocal dilution buffer**

Store at 4 °C. Use within one month.

Reagent	Final concentration	Quantity or Volume
---------	---------------------	--------------------

PBS (Recipe 3)	1×	10 mL
BSA	1% (w/v)	0.1 g
Triton X-100	0.25% (v/v)	25 µL
Total	n/a	10 mL

## Laboratory supplies

1. Vacuette<sup>®</sup> tube 4 mL CAT serum clot activator blood collection tubes (Greiner, Bio-one, catalog number: 454204)
2. Vacuette<sup>®</sup> tube 9 mL NH sodium heparin blood collection tubes (Greiner, Bio-one, catalog number: 455051)
3. Microcentrifuge tubes, 0.5 mL, 1.5 mL, 2 mL (Greiner, Bio-one, catalog numbers: 667201, 616201, 623201)
4. 1.5 mL microcentrifuge tubes, RNase free (SARSTEDT, catalog number: 72.706.700)
5. PCR tubes and caps (VWR, catalog number: 20170-010)
6. 50 mL conical centrifuge tubes (Greiner Bio-One, catalog number: 227261)
7. Serological pipettes, 2 mL, 5 mL, 10 mL, 50 mL (Greiner Bio-One, catalog numbers: 710180, 606180, 607180, 760180, 768160)
8. Screw cap microtubes (SARSTEDT, catalog number: 72.694.305)
9. Microlance<sup>™</sup> 3 needles, 21G (BD, catalog number: 304432)
10. 2.5 mL syringe (Greiner Bio-One, catalog number: SYR2)
11. 8.5 cm Petri dishes (Greiner, Bio-one, catalog number: 633181)
12. Pipette tips, filtered, 20 µL, 200 µL, 1,000 µL (Mettler Toledo, Ranin, catalog numbers: 17005861, 17005863, 30389211)
13. Gene Pulser electroporation cuvettes, 0.1 cm gap (Bio-Rad, catalog number: 1652083)
14. Polystyrene semi micro cuvettes (Fisherbrand, catalog number: FB55147)
15. 125 mL baffled flasks (VWR, catalog number: 214-0458)
16. 250 mL baffled flasks (VWR, catalog number: 214-0460)
17. 1 L baffled flasks (VWR, catalog number: 214-0464)
18. Micro quartz absorption cuvettes (Merck, Hellma<sup>®</sup>, catalog number: Z600210)
19. Square bioassay dish (Nunc, catalog number: 240845)
20. Reagent reservoirs, non-sterile (VWR, catalog number: 613-1176)
21. Reagent reservoirs, sterile (Axygen, catalog number: RES-V-25-S)
22. Deep-well 96-well 2 mL plate (Fisherbrand, catalog number: 11391555)
23. Deep-well 96-well 1 mL plate (Greiner, Bio-one, catalog number: 780215)
24. Deep-well 96-well 0.5 mL plate (Greiner, Bio-one, catalog number: 786261)
25. Cell culture adhesive seal (Azenta Life Sciences, catalog number: 4ti-0517)
26. Adhesive film for ELISA and general incubation (VWR, catalog number: 60941-062)
27. PCR foil seal (Azenta Life Sciences, catalog number: 4ti-0550)
28. Rigid semi-skirted 96-well PCR plate (Thermo Scientific, catalog number: AB-0990)
29. Framestar<sup>®</sup> 96-well skirted PCR plate (Azenta Life Sciences, catalog number: 4ti-0960)
30. 96-well ELISA microplate (Greiner Bio-One, catalog number: 655061)
31. Octet SA biosensors (Sartorius, catalog number: 18-5019)
32. 96-well black plate, polypropylene (Greiner Bio-One, catalog number: 655209)
33. 96-well microtitre plate (Greiner Bio-One, catalog number: 65101)
34. 250 mL PPCO centrifuge tube (Nalgene, catalog number: 3141-0250)
35. 25 mL high speed PPCO centrifuge tube (Nalgene, catalog number: 3119-0010)
36. NuPage 4%–12% bis tris precast gel (Invitrogen, catalog number: NP0323BOX). A suitable alternative is homemade 12% running and 4% stacking SDS-PAGE gel
37. 0.22 µm syringe filters (Agilent, catalog number: 5190-5116)
38. µ-Slide 8-well high polymer bottom chambered coverslip (ibidi, catalog number: 80807)
39. Microscope slides with ground edges (Fisherbrand, catalog number: 11572203)
40. Microscope coverslip, 12 mm diameter, no. 1.5 (Scientific Laboratory supplies, catalog number: MIC3334)
41. 6-well cell culture dish (Greiner Bio-One, catalog number: 657160)

42. Dumont tweezer style 5 (Electron Microscopy Sciences, catalog number: 0203-5-PS)

## Equipment

1. Pipetboy acu 2 pipette controller (Ingetra, catalog number: I155017)
2. Single-channel pipettes starter kit (Mettler Toledo, Ranin, catalog number: 30386738)
3. Multichannel pipettes 1–10  $\mu$ L, 20–200  $\mu$ L, 100–1,200  $\mu$ L (Mettler Toledo, Ranin, catalog numbers: 17013802, 17013805, 17014496)
4. Microplate washer (Thermo Fisher, Wellwash™, catalog number: 5165000). If there is no well washer available, the wells can be washed manually using a multichannel pipette or using a wash bottle filled with PBST
5. CLARIOstar plus plate reader (BMG Labtech, catalog number: 430-501S-FL)
6. Class 2 microbiological safety cabinet (referred to as a safety cabinet and known as a laminar flow) (Contained Air solutions, catalog number: BioMAT 2). If no safety cabinet is available, working by a Bunsen burner is a suitable alternative. Working by an open flame with flammable EtOH should be done with caution
7. PCR workstation cabinet (Bigneat, catalog number: MW520-20)
8. NanoDrop™ One/One<sup>C</sup> Microvolume UV-Vis spectrophotometer (Thermo Scientific, catalog number: ND-ONE-W)
9. DynaMag™ 96-side magnet (Invitrogen, catalog number: 12331D)
10. PCR thermal cycler (Applied Biosystems, Veriti™, catalog number: 4375305)
11. Eporator (Eppendorf, catalog number: 4309000027)
12. ChemiDoc™ imaging system (Bio-Rad, catalog number: 12003154)
13. Cell density meter (Fisherbrand, catalog number: A0)
14. Multi Sub Electrophoresis System (Flowgen, catalog number: FMMS10)
15. SDS-PAGE equipment, mini gel tank (Invitrogen, catalog number: A25977)
16. PowerPac basic power supply (Bio-Rad, catalog number: 1645050EDU)
17. JB Nova water bath (Grant, catalog number: JBN5)
18. Dry heating block (Grant, catalog number: QBD2)
19. Dual LED Blue/White Light Transilluminator (Thermo Fisher, catalog number: LB0100)
20. Safe Imager™ viewing glasses (Thermo Fisher, catalog number: S37103)
21. HulaMixer™ sample mixer (Invitrogen, catalog number: 15920D). Agitation using this sample mixer is performed using the following settings: orbital = 5, rpm = 1
22. DynaMag™-2 magnet (Invitrogen, catalog number: 12321D)
23. LSE™ digital microplate shaker (Corning, catalog number: 4782-4). Agitation using the microplate shaker is performed at 500 rpm
24. Multifuge X4 Pro-MD (Thermo Scientific, catalog number: 75009500). In all centrifugation steps, maximum acceleration and deceleration rates are used unless otherwise specified
25. Fresco™ 21 microcentrifuge (Thermo Scientific, catalog number: 75002555)
26. Orbital shaking incubator (Shel Lab, catalog number: SI6/SI6R) (for 600 rpm agitation steps)
27. 44R incubator shaker (New Brunswick, Innova®, catalog number: M1282-0006) (for 200 and 250 rpm agitation steps)
28. Octet® R8 (Sartorius, catalog number: 30-0518)
29. SP8 Lightning confocal microscope (Leica Microsystems Ltd.)
30. Heracell™ VIOS 160i CO<sub>2</sub> Incubator (Thermo Scientific, catalog number: 51033559)
31. Stuart™ SSM3 Gyratory Rocker (Cole-Parmer, catalog number: 51900-26)

## Software and datasets

1. SnapGene (version 7.1.0, 2023, released 28 November 2023); requires a license
2. Octet® Analysis Studio Software (version 12.2.2.26); requires a license

3. LasX software for confocal microscope (version 5.2.1); requires a license

## Procedure

### A. Llama immunisation

Antibodies are raised in a llama by intramuscular immunisation with up to eight different proteins in parallel. The identity of the proteins depends on the intended application. We have raised nanobodies to integral membrane proteins (e.g., PEPT2 [4]), cell surface glycoproteins (e.g., GPC3 [7]) and viral antigens (e.g., SARS-CoV-2 spike protein [6]).

1. Store protein antigens at -80 °C in  $3 \times 200 \mu\text{g}$  aliquots, preferably  $200 \mu\text{L}$  of  $1 \text{ mg/mL}$  per aliquot, in PBS or a buffer that is optimal for the protein (e.g., if in a detergent-solubilised membrane protein, then the required amount of appropriate detergent would be included in the buffer).
2. Before immunisation, collect  $5 \text{ mL}$  of blood in blood tubes without anticoagulant to prepare a sample of pre-immune serum. Preparation of the pre-immune serum sample is detailed in step C2.
3. Thaw an aliquot of each of the proteins that are to be included in the immunisation. Gently mix each protein with an equal volume of Gerbu adjuvant and inject subcutaneously (maximum  $2 \text{ mL}$  per site) in the neck base/shoulder of the llama.
4. The llama is immunised on day 0 followed by two boosts of antigen on days 28 and 56. Ten days following the final boost,  $170 \text{ mL}$  of blood from the jugular vein is collected into heparinised tubes (to prevent coagulation) for isolation of polymorphonuclear cells (PMNCs). In addition, collect  $5 \text{ mL}$  of blood in blood tubes without anticoagulant to prepare a sample of post-immune serum. Details for how to prepare the post-immune serum sample are described in step C2.

### B. Seroconversion ELISA

A seroconversion ELISA is carried out to confirm if the immunisation of the llama was successful in generating a response to the injected antigen(s). A strong ELISA signal is strongly suggestive that binding nanobodies can be isolated from the created library by phage display.

1. Coat a 96-well ELISA microplate with  $100 \mu\text{L/well}$  of  $10 \mu\text{g/mL}$  neutrAvidin™ diluted in PBS and incubate overnight at  $4^\circ\text{C}$ .
2. Using a plate washer, wash the plate five times with  $300 \mu\text{L/well}$  of PBST.
3. Add  $100 \mu\text{L/well}$  of  $50 \text{ nM}$  biotinylated target protein diluted in PBS and incubate for  $1 \text{ h}$  at RT on a microplate shaker. See **General note 1**.
4. Using a plate washer, wash the plate five times with  $300 \mu\text{L/well}$  of PBST.
5. Add  $250 \mu\text{L/well}$  of blocking solution and incubate for  $1 \text{ h}$  at RT on a microplate shaker.
6. Using a plate washer, wash the plate five times with  $300 \mu\text{L/well}$  of PBST.
7. Add  $110 \mu\text{L}$  of either pre- or post-immunisation serum diluted 1:10 in blocking solution to the first well of the 96-well plate. Prepare a serial dilution of the sera in the plate by adding  $100 \mu\text{L}$  of blocking solution to the second well and add  $10 \mu\text{L}$  of the solution from the first well. Repeat for a third well so that the sera are diluted 1 in 10, 100, and 1,000; include a no-serum control. Incubate for  $1 \text{ h}$  at RT on a microplate shaker.
8. Using a plate washer, wash the plate five times with  $300 \mu\text{L/well}$  of PBST.
9. Add  $100 \mu\text{L/well}$  of anti-llama-HRP diluted 1:2,500 in 0.1% BSA-PBS and incubate for  $1 \text{ h}$  at RT on a microplate shaker.
10. Using a plate washer, wash the plate five times with  $300 \mu\text{L/well}$  of PBST.
11. Add  $100 \mu\text{L/well}$  of ABTS substrate, prepared by mixing solution A and solution B in a 1:1 ratio. Protect from light and measure the absorbance at  $405 \text{ nm}$  after  $15 \text{ min}$ . See **General note 2**.

### C. VHH library preparation

Below (Table 1) is a suggested timetable of experiments required for the preparation of a VHH library. Instances where the protocol can be interrupted/held are indicated by “pause point.”

**Table 1. Experimental timetable for the preparation of a VHH library**

	Task 1	Task 2	Task 3
Day 1	Isolation of PMNCs (step C1)	Isolation of immune serum (step C2)	RNA extraction (steps C3a–g) <b>Pause point</b>
Day 2	RNA extraction (steps C3h–i)	Reverse transcription and PCR amplification of VHH1 (step C4)	Production of VHH2 (step C5) <b>Pause point</b>
Day 3	Digestion of pADL23c (step C6) <b>Pause point</b>		
Day 4	Test of library size (steps C7a–g)		
Day 5	Test of library size (steps C7h–n) <b>Pause point</b>		
Day 6	Scaled library preparation (steps C8a–i)		
Day 7	Scaled library preparation (steps C8j) <b>Pause point</b>		

1. Isolation of PMNCs

The below instructions are based on receiving 170 mL of llama blood in 17 × 10 mL blood tubes without anticoagulant. Isolation of the PMNCs should be done on the day of the blood draw. All work is performed in a safety cabinet after spraying with 5% Chemgene followed by 70% EtOH.

**Caution:** Decontaminate all blood-contaminated consumables and liquid with 2% Virkon, followed by autoclaving.


- Transfer blood from 2 × 10 mL blood tubes into 1 × 50 mL conical centrifuge tube using a 10 mL serological pipette.
- Add 10 mL of PBS to each of the empty blood vials and transfer into the conical centrifuge tube with the 20 mL of blood (step C1a) to give a total volume of 40 mL of diluted blood per 50 mL conical centrifuge tube, or a total of 340 mL of diluted blood. Invert tubes gently to ensure a homogenous mixture.
- In fresh 23 × 50 mL conical centrifuge tubes, add 15 mL of Histopaque®-1077 and then gently layer 15 mL of the diluted blood on top. See **General note 3**. Hold the tube at a 45° angle from the horizontal while adding the blood. See **General note 4**.

**Critical:** Avoid agitating the tubes at this point to maintain the density gradient that is forming.

- Centrifuge the conical centrifuge tubes at 800× g for 20 min at 18 °C using an acceleration = 1 and deceleration = 0.

**Critical:** The acceleration and deceleration rates are very important, as they will create and maintain the desired gradient. The results of the density gradient and localisation of the PMNC layer is shown in Figure 1.

- From the 23 prepared blood–Histopaque®-1077 tubes, collect the PMNC layer at the plasma–Histopaque®-1077 interface (red square in Figure 1). This is achieved by using a 25 mL serological pipette and placing it just above the interface, drawing the liquid up whilst moving across the surface interface. Any Histopaque®-1077 or serum that is drawn up with the PMNCs will be washed away in the next steps.

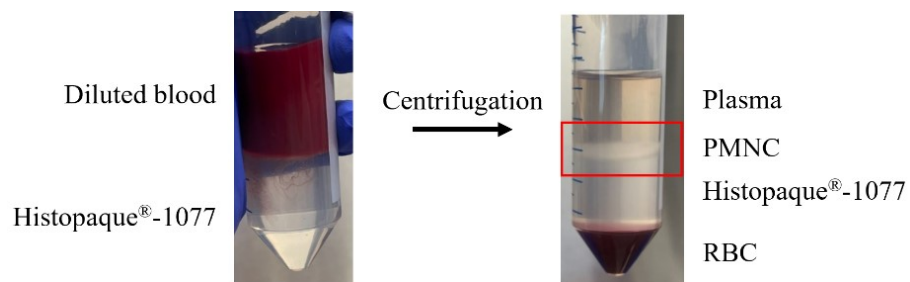

**Figure 1. Isolation of polymorphonuclear cells (PMNCs) from llama blood using a density gradient.** Diluted llama blood overlaid on Histopaque®-1007 and the resulting layer of PMNCs and red blood cells (RBCs) after centrifugation.

- f. Transfer 2× conical centrifuge tubes' worth of isolated PMNCs into a single 50 mL conical centrifuge tube. Typically, 15–20 mL of PMNCs is retrieved from each conical centrifuge tube. If more is collected, ensure that 20 mL of isolated PMNCs is added to fresh 50 mL conical centrifuge tubes. There should be between 23 and 25 × 50 mL conical centrifuge tubes at this point.
- g. To each conical centrifuge tube containing isolated PMNCs, add PBS to a total volume of 40 mL. One can just pour from the bottle of PBS, as long as it does not touch the conical centrifuge tube. Gently invert the conical centrifuge tubes several times and centrifuge at 250× g for 10 min at 18 °C.
- h. Pour off the supernatant and gently resuspend the pellet from each conical centrifuge tube in 5 mL PBS using a new 5 mL serological pipette per pellet. Centrifuge the cells at 100× g for 5 min at 18 °C.

*Note: At this point, start using RNaseZap™ to wipe the work surface, gloves, and any items that are brought into the safety cabinet.*

- i. Use a 5 mL serological pipette to remove the supernatant from each 50 mL conical centrifuge tube. Gently resuspend the pellets from three conical centrifuge tubes in a total of 1 mL of PBS into a single 50 mL conical centrifuge tube using a 2 mL serological pipette. Use the final amount of liquid in this serological pipette and dispense it into a 1.5 mL microcentrifuge tube. There should be six 50 mL conical centrifuge tubes and six 1.5 mL microcentrifuge tubes at this point.
- j. Remove the six 1.5 mL microcentrifuge tubes from the safety cabinet. In the lid of each of the 1.5 mL microcentrifuge tubes, add 20 µL of cells and 20 µL of trypan blue and mix gently. Determine cell number and viability using an automated cell counter. Expect approximately  $1 \times 10^7$  cells at 98%–100% viability per sample taken from each conical centrifuge tube. Thus, an overall yield of  $6 \times 10^7$  cells should be expected. See **General note 5**.

*Note: The following steps can be performed at a lab bench or in a PCR cabinet after being sprayed with RNaseZap™ and using RNase-free tips and 1.5 mL microcentrifuge tubes.*

- k. Transfer the cell suspension from each 50 mL conical centrifuge tube (step C1i) into their own 1.5 mL microcentrifuge tube. There should be six 1.5 mL microcentrifuge tubes in total. Pellet the cells at 100× g for 5 min at 4 °C. Discard the supernatant. The RNA can now be isolated from the pellets (step C3a).

**Critical:** Perform RNA extraction (step C3) immediately.

2. Isolation of pre- and post-immune serum  
Isolation of serum must be done on the day of the blood draw. These steps can be done during the isolation of PMNCs at a lab bench.
  - a. Incubate the coagulated blood at RT for 30 min.
  - b. Transfer the serum and the blood clot into a 50 mL conical centrifuge tube. This is achieved by using a 2 mL serological pipette to stab holes in the blood clot and using the stripette to guide the clot and serum into a 50 mL conical centrifuge tube. Centrifuge the conical centrifuge tube at 1,000× g for 10 min at 4 °C.
  - c. Transfer the serum into a 15 mL conical centrifuge tube. A second centrifugation at 1,000× g for 10



min at 4 °C may be required to obtain clear serum. Prepare 500 µL aliquots of serum in screw cap microtubes, which are then stored at -20 °C.

### 3. RNA extraction

*Note: All work should be done at a lab bench or in a PCR cabinet after spraying with RNaseZap™. Use RNase-free consumables and spray gloves with RNaseZap™. See General note 6.*

- To each cell pellet in each of the six 1.5 mL microcentrifuge tubes (step C1k), add 1 mL of cold TRIzol™ reagent. Using a 21 G needle and 2.5 mL syringe, homogenise the contents of each 1.5 mL microcentrifuge tube by aspirating and dispensing ten times through the needle.

**Caution:** Use caution when working with needles.

- Incubate the resulting lysate on ice for 5 min to permit complete dissociation of the nucleoproteins complex and then add 0.2 mL of chloroform to each 1.5 mL microcentrifuge tube.
  - Vortex each microcentrifuge tube vigorously for approximately 10 s. A homogenous pink milky solution should occur. Incubate on ice for 3 min and then centrifuge at 12,000× g for 15 min at 4 °C. The mixture separates into a lower red phenol-chloroform layer, an interphase layer, and a colourless upper aqueous phase. If the upper layer is still cloudy or has not formed properly, vortex and centrifuge again.
  - Transfer approximately 500 µL of the upper aqueous phase containing the RNA to a new 1.5 mL microcentrifuge tube using a 200 µL pipette. Some of the upper aqueous layer can be left behind to ensure that none of the precipitate from the interphase is included. Add 0.5 mL of isopropanol, gently invert each tube several times, and incubate on ice for 10 min.
  - Centrifuge at 12,000× g for 10 min at 4 °C to pellet the RNA. Keep the lid of the 1.5 mL microcentrifuge tube facing outwards in the rotor to aid in the visualisation of the white gel-like pellet after centrifugation. Remove the approximately 900 µL of the isopropanol containing supernatant using a 200 µL pipette to ensure that the pellet is not disturbed or accidentally aspirated with the waste.
  - Add 1 mL of 75% ethanol to the 1.5 mL microcentrifuge tube, vortex briefly, and centrifuge at 7,500× g for 5 min at 4 °C. The RNA pellet will look a little whiter and slightly bigger than what it did at the end of the previous step. Remove the approximately 900 µL of the supernatant using a 200 µL pipette to ensure that the pellet is not disturbed or accidentally aspirated.
  - Add 1 mL of 75% EtOH and store the purified RNA at -20 °C. See General note 7.
- Pause point.**
- Remove one microcentrifuge tube of RNA from the -20 °C and centrifuge at 7,500× g for 5 min at 4 °C. Discard the EtOH using a 200 µL and then a 10 µL pipette to remove all EtOH without accidentally aspirating the pellet. Allow to dry at RT for 10 min.
  - Dissolve the RNA pellet in 30 µL of nuclease-free water and measure the A<sub>260</sub> using the Nanodrop. Expect a concentration between 500 and 800 ng/µL per 1.5 mL tube with values of approximately 2.0 for both A<sub>260/280</sub> and A<sub>260/230</sub>.

**Critical:** Use RNase-free tips and pipette from the PCR cabinet to measure the RNA concentration to avoid the introduction of RNases from the communal tips and pipette used at the Nanodrop. Once the RNA has been isolated, proceed to step C4 as soon as possible to produce cDNA, which is more stable than RNA.

### 4. Reverse transcription and PCR amplification of VHH1

*Note: Work at a lab bench or in a PCR cabinet after spraying with RNaseZap™. Use RNase-free consumables and spray gloves with RNaseZap™. See General note 6.*

- Prepare the following reaction in a PCR tube (1 × 50 µL) (Table 2).

**Table 2. PCR reaction composition**

Reagent	Final concentration	Volume
Nuclease-free water	n/a	9.5 µL
10 µM CALL 001 primer	0.5 µM	2.5 µL
10 µM CALL 002 primer	0.5 µM	2.5 µL
2× Platinum™ SuperFi™ RT-PCR Master Mix	1×	25 µL
1 µg Template RNA (step C3i)	200 ng	10 µL

SuperScript™ IV RT Mix		0.5 µL
Total	n/a	50 µL

- b. Briefly vortex and centrifuge the reaction tube and place on ice.  
c. Set up the thermal cycler as below and perform the PCR (Table 3).

Table 3. PCR cycling conditions

Step	Temperature (°C)	Duration	Number of cycles
Reverse transcription	55	10 min	
Reverse transcription inactivation/initial denaturation	98	2 min	1
Denaturation	98	10 s	
Annealing	66	10 s	30
Extension	72	30 s	
Final extension	72	5 min	1
Hold	12	Infinite hold	-

*Note: There is no need to maintain RNase-free conditions from this point onwards. Work at a lab bench.*

- d. Add 10 µL of 6× DNA gel loading buffer to the reaction tube and run alongside a lane of 1× GeneRuler 1 kb DNA ladder on a 0.7% agarose gel containing 1× SYBR™ safe DNA gel stain. Electrophorese in 1× TBE buffer at 80 V for 40 min.
- e. Visualise the amplified 700 bp fragment (VHH-CH2) using a LED Blue transilluminator, excise from the gel using a scalpel knife, and gel extract using the Nucleospin Gel and PCR clean-up kit according to the manufacturer’s protocol. Elute in a final volume of 50 µL of NE buffer.  
**Caution:** Use Safe Imager™ viewing glasses when visualising the DNA during excision of the band from the agarose gel and use caution when using scalpel knives.
- f. Using the Purelink PCR purification kit, add 4× volume of B2 buffer (200 µL) to the gel-extracted product and purify using one Purelink column from the Purelink PCR purification kit eluting in 25 µL E1 buffer. Measure the A<sub>280</sub> and expect a concentration of approximately 40–50 ng/µL.
- g. Mix 6× DNA gel loading buffer with 120 ng of VHH1 (step C4f) to a final concentration of 1× and run alongside a well of 1× HyperLadder™ 1 kb DNA ladder on a 1% agarose gel containing 1× SYBR™ safe DNA gel stain to determine the purity. Electrophorese in 1× TBE buffer at 80 V for 40 min. Expect a result as shown in Figure 2.

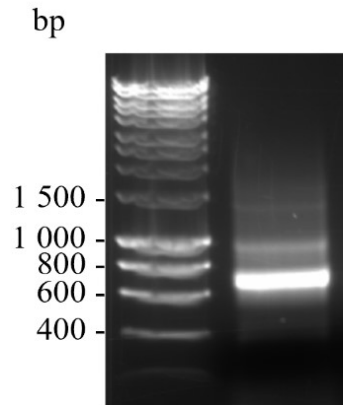


Figure 2. Sample of amplified VHH1 of approximately 700 bp visualised on a 1% agarose gel

containing 1× SYBR™ safe DNA gel stain. The sizes in base pairs of a DNA ladder run in parallel are shown.

5. Production of VHH2
  - a. Prepare the following PCR reaction in PCR tubes (6 × 50 µL reactions) (Table 4).

**Table 4. PCR reaction composition**

Reagent	Final concentration	Volume
Nuclease-free water	n/a	19 µL
10 µM VHH For2 primer	0.5 µM	2.5 µL
10 µM VHH Rev primer	0.5 µM	2.5 µL
2× Phusion flash PCR master mix	1×	25 µL
5 ng/µL VHH1 (step C4f)	1 ng	10 µL
Total	n/a	50 µL


- b. Briefly vortex and centrifuge the reaction tubes.
- c. Set up the thermal cycler as below and perform the PCR (Table 5).

**Table 5. PCR cycling conditions**

Step	Temperature (°C)	Duration	Number of cycles
Denaturation	98	30 s	1
Annealing	98	10 s	
Extension	55	30 s	30
Final extension	72	20 s	
Hold	72	5 min	1
Denaturation	12	Infinite hold	-

- d. Combine all six reactions into a single 2 mL microcentrifuge tube. Using the Purelink PCR purification kit, add 4× volume of B2 buffer (1.2 mL) to the amplified PCR product and purify between two Purelink columns, eluting each with 25 µL of E1 buffer. Combine both eluates into one 1.5 mL microcentrifuge tube, measure the A<sub>280</sub>, and expect a concentration between 200 and 300 ng/µL.
- e. Mix 6× DNA gel loading buffer with 120 ng of VHH2 (step C5d) to a final concentration of 1× and run alongside a well of 1× Hyperladder™ 1 kb DNA ladder on a 1% agarose gel containing 1× SYBR™ safe DNA gel stain to determine the purity. Electrophoresis in 1× TBE buffer at 80 V for 40 min. Expect a result as shown in Figure 3.

**Pause point.**

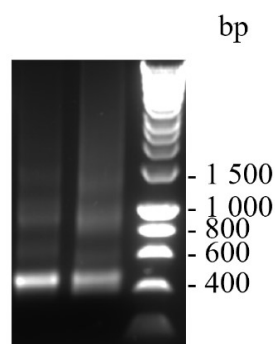

**Figure 3. Sample of amplified VHH2 of approximately 400 bp visualised on a 1% agarose gel containing 1× SYBR™ safe DNA gel stain. The sizes in base pairs of a DNA ladder run in parallel**

are shown.


6. Digestion of pADL23c

- a. Prepare the following restriction digestion reaction in PCR tubes ( $10 \times 100 \mu\text{L}$  reactions) (Table 6).

**Table 6. Restriction digestion reaction composition**

Reagent	Final concentration	Volume
Nuclease-free water	n/a	83 $\mu\text{L}$
rCutSmart <sup>TM</sup> buffer	1×	10 $\mu\text{L}$
SfiI	5 units	2.5 $\mu\text{L}$
400 ng/ $\mu\text{L}$ vector DNA	1.8 $\mu\text{g}$	4.5 $\mu\text{L}$
Total	n/a	100 $\mu\text{L}$


- b. Briefly vortex and centrifuge the reaction tubes.  
c. Incubate at 50 °C for 2 h.  
d. Add 20  $\mu\text{L}$  of 6× DNA gel loading buffer to each reaction tube and run alongside a well of 1× GeneRuler 1 kb DNA ladder on a 0.7% agarose gel containing 1× SYBR<sup>TM</sup> Safe DNA gel stain. Electrophorese in 1× TBE buffer at 80 V for 40 min. Expect a result as shown in Figure 4.

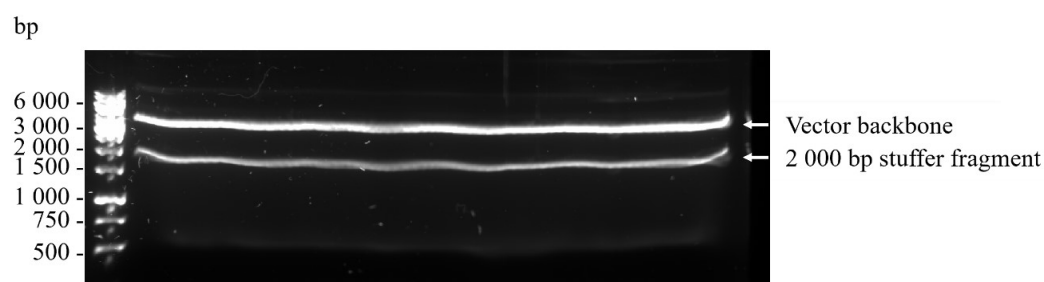

**Figure 4. Preparative 0.7% agarose gel containing 1× SYBR<sup>TM</sup> safe DNA gel stain of the vector backbone (upper band) and stuffer fragment (lower band) after SfiI digestion of pADL23c phagemid vector prior to gel extraction.** The sizes in base pairs of a DNA ladder run in parallel are shown.

- e. Visualise the digested, higher molecular weight, approximately 3,000 bp vector backbone using a LED Blue transilluminator, excise from the gel using a scalpel knife, and gel extract using 10 spin columns from the Nucleospin Gel and PCR clean-up kit according to the manufacturer's protocol. Elute each column using 50  $\mu\text{L}$  of NE buffer and pool all the eluates into a single 1.5 mL microcentrifuge tube.  
**Caution:** Use Safe Imager<sup>TM</sup> viewing glasses when visualising the DNA during excision of the band from the agarose gel and use caution when using a scalpel knife.  
f. Using the Purelink PCR purification kit, add 4× volume of B3 buffer (2 mL) to the gel extracted product and purify using one Purelink column from the Purelink PCR purification kit eluting in 50  $\mu\text{L}$  E1 buffer. Measure the  $A_{280}$  and expect a concentration between 50 and 80 ng/ $\mu\text{L}$ .  
g. Mix 6× DNA gel loading buffer with 120 ng of VHH2 (step C5d) and digest pADL23c backbone (step C6f) to a final concentration of 1× and run alongside a well of 1× HyperLadder<sup>TM</sup> 1 kb DNA ladder on a 1% agarose gel containing 1× SYBR<sup>TM</sup> Safe DNA gel stain to determine the purity. Electrophorese in 1× TBE buffer at 80 V for 40 min. Expect a result as shown in Figure 5.

**Pause point.**

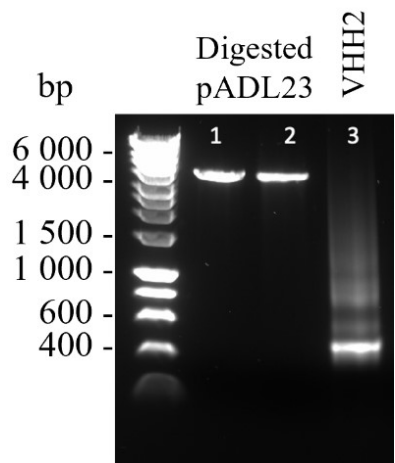

**Figure 5.** Samples of linearised and purified pADL23c (lanes 1 and 2) and purified VHH2 (lane 3) visualised on a 1% agarose gel containing 1× SYBR™ safe DNA gel stain. The sizes in base pairs of a DNA ladder run in parallel are shown.

7. Test of library size
- a. Prepare the following in-fusion reaction in a PCR tube (1 × 10 µL) (Table 7).

**Table 7. PCR reaction composition**

Reagent	Final concentration	Volume
Nuclease-free water	n/a	5 µL
5× CEII buffer	1×	2 µL
100 ng/µL VHH2 insert (step C5d)	100 ng	1 µL
20 ng/µL vector DNA (step C6f)	20 ng	1 µL
Exnase II		1 µL
Total	n/a	10 µL

- b. Briefly vortex and centrifuge the reaction tube.
- c. Incubate at 42 °C for 30 min using a thermal cycler. Add 40 µL of TE buffer to terminate the reaction.
- d. Using the Purelink PCR purification kit, add 4× volume of B3 buffer (200 µL) to the in-fusion product and purify using one Purelink column eluting in 20 µL nuclease-free water.

*Note: Work in a safety cabinet.*

- e. Gently mix 5 µL of purified in-fusion reaction with 30 µL of competent TG1 in a 1.5 mL microcentrifuge tube and transfer to a chilled 1 mL electroporation cuvette.
- f. Electroporate at 1.7 kV and immediately add 400 µL of recovery medium to the cuvette. Transfer contents to a 2 mL microcentrifuge tube and incubate for 1 h at 37 °C agitating at 600 rpm. See **General note 8**.
- g. Using 50 µL of culture, prepare a 1:10, 1:100, and 1:1,000 dilution in 2× YT. Spread 100 µL on 1% agar LB plates containing 100 µg/mL ampicillin until dry. Incubate overnight at 37 °C.
- h. Estimate the library size in terms of colony forming units per millilitre (CFU/mL) using the following formula: number colonies × 10 (to get to mL) × 0.4 (volume) × dilution. See **General note 9**.
- i. Based on the estimated small-scale library size, set up the required number of reactions considering that 5 µL of in-fusion reaction yielded × size library. Calculate how many microlitres of in-fusion reactions would be needed to make a library with a size of 1 × 10<sup>6</sup>. See **General note 10**.
- j. To validate the number of full-length VHH clones in the small-scale library, a colony PCR to amplify the ~500 bp VHH gene in 48 randomly selected clones from a plate in step C7g is performed. Prepare the following PCR reaction master mixture in a 2 mL microcentrifuge tube (Table 8).

Table 8. PCR reaction composition

Reagent	Final concentration	Volume for single reaction	Volume for master mixture
Nuclease-free water	n/a	20.75 µL	1037.5 µL
10 µM PhD seq Fwd primer	0.2 µM	0.5 µL	25 µL
10 µM PhD seq Rev primer	0.2 µM	0.5 µL	25 µL
10× Taq buffer	1×	2.5 µL	125 µL
25 mM dNTP	0.25 mM	0.25 µL	12.5 µL
Taq polymerase	1 unit	0.5 µL	25 µL
Total	n/a	25 µL	1250 µL

- k. Briefly vortex and centrifuge the reaction tube.
- l. Pipette 25 µL of the master mixture into 48 wells of a 96-well rigid semi skirted PCR plate. Pick 48 colonies with 48 × 10 µL tips and place into each filled well. Remove the tips using a multichannel pipette. Cover the PCR plate with PCR foil seal.
- m. Set up the thermal cycler as below and perform the PCR (Table 9).

Table 9. PCR cycling conditions

Step	Temperature (°C)	Duration	Number of cycles
Denaturation	95	7 min s	1
Annealing	95	15 s	
Extension	55	30 s	35
Final extension	68	1 min 40 s	
Hold	68	5 min	1
Denaturation	12	Infinite hold	-

- n. Add 20 µL of 2× DNA gel loading buffer to each well of the PCR plate and load 10 µL to each well alongside a well of 1× GeneRuler 1 kb DNA ladder on a 1% agarose gel containing 1× SYBR™ safe DNA gel stain. Electrophoresis in 1× TBE buffer at 80 V for 40 min. In Figure 6, 46 out of the 48 clones that were subjected to colony PCR possessed the ~500 bp band corresponding to the amplified VHH gene. The two clones that were not positive are indicated by an asterisk.

Pause point.

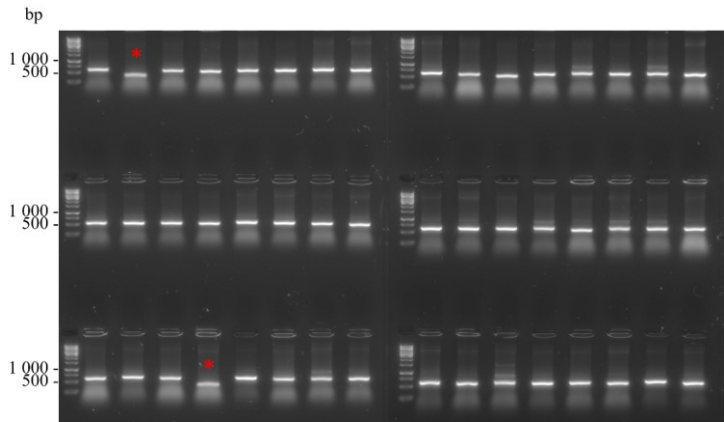


Figure 6. Colony PCR of 48 randomly selected clones from a small-scale library preparation

8. Scaled library preparation
  - a. Prepare the following reaction in PCR tubes ( $20 \times 20 \mu\text{L}$ ) (Table 10).


**Table 10. In-Fusion reaction composition**

Reagent	Final concentration	Volume
Nuclease-free water	n/a	9 $\mu\text{L}$
5 $\times$ CEII buffer	1 $\times$	4 $\mu\text{L}$
50 ng/ $\mu\text{L}$ VHH2 insert (step C5c)	100 ng	2 $\mu\text{L}$
5 ng/ $\mu\text{L}$ vector DNA (step C6f)	20 ng	4 $\mu\text{L}$
Exnase II		1 $\mu\text{L}$
Total	n/a	20 $\mu\text{L}$

- b. Briefly vortex and centrifuge the reaction tube.
- c. Incubate at 42 °C for 30 min in a thermal cycler.
- d. Combine reactions into a single 2 mL microcentrifuge tube. Using the Purelink PCR purification kit, add 4 $\times$  volume of B3 buffer (1.6 mL) to the in-fusion product and purify using one Purelink column eluting in 50  $\mu\text{L}$  of nuclease-free water. Measure the  $A_{280}$  and expect a concentration of between 35 and 50 ng/ $\mu\text{L}$ .

*Note: Work in a safety cabinet.*

- e. TG1 cells are transformed by electroporation in separate reactions to create the scaled library. Gently mix 6  $\mu\text{L}$  of purified in-fusion reaction with 30  $\mu\text{L}$  of competent TG1 in a 1.5 mL microcentrifuge tube and transfer to a chilled 1 mL electroporation cuvette. Electroporate at 1.7 kV, immediately add 960  $\mu\text{L}$  of recovery medium to the cuvette, and transfer to a 50 mL conical centrifuge tube. See **General note 8**.
- f. Repeat this a total of eight times with eight separate in-fusion-TG1 mixtures.
- g. Combine two cuvettes worth of electroporated cells into a single 50 mL conical centrifuge tube and incubate at 37 °C for 1 h shaking at 200 rpm. You should have four conical centrifuge tubes in total.
- h. Use 50  $\mu\text{L}$  from one of the conical centrifuge tubes and prepare a 1:10, 1:100, and 1:1,000 dilution in 2 $\times$  YT. Spread 100  $\mu\text{L}$  of each dilution on 1% agar LB plates containing 100  $\mu\text{g}/\text{mL}$  ampicillin until dry. Incubate overnight at 37 °C. Use these plates to estimate library size ( $\text{CFU}/\text{mL}$ ) = number colonies  $\times 10$  (to get to mL)  $\times 1.92$  (volume)  $\times$  dilution. See **General note 9**.
- i. Plate the entire culture of each of the four conical centrifuge tubes onto separate bioassay dishes filled with 2% agar LB plates containing 100  $\mu\text{g}/\text{mL}$  ampicillin. There should be four plates in total. Spread the entire volume until dry and incubate overnight at 37 °C. Expect a density similar to that shown in Figure 7, which is equivalent to a  $1 \times 10^6$  CFU/mL sized library. See **General note 11**.

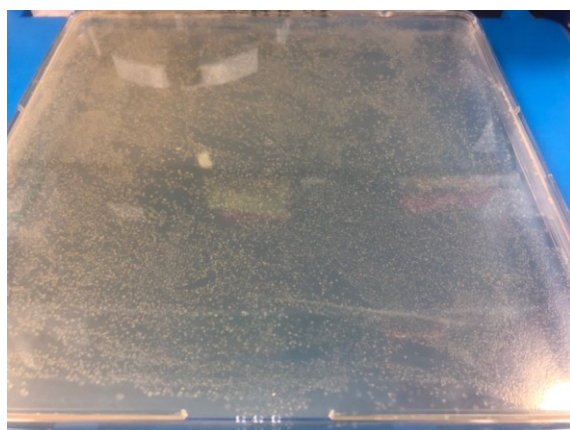

**Figure 7. Expected density on a bioassay plate for a  $1 \times 10^6$  CFU/mL library**



- j. Resuspend the lawn of colonies of each bioassay plate in 6 mL of 2× YT (containing a final concentration of 25% glycerol) and transfer to a single 50 mL conical centrifuge tube. Mix well by serological pipette and prepare 1 mL aliquots in screw cap microtubes, which are stored at -80 °C. This is the VHH library stock, which is used for the first round of panning.

**Pause point.**

## D. Panning

Below is a suggested timetable for panning, hit identification by ELISA, and small-scale expression to maximise time in the lab and to allow space for overlap if multiple panning campaigns are being performed for multiple proteins by multiple people using shared resources (Table 11).

1. Prepare log-phase TG1

*Note: Work in a safety cabinet.*

- a. Streak TG1 on a 1% agar LB plate without antibiotics and incubate overnight at 37 °C.
- b. Inoculate 20 mL of 2× YT with a single TG1 colony in a 100 mL disposable baffled Erlenmeyer flask and incubate overnight at 37 °C shaking at 200 rpm.
- c. Mix 20 mL of 50% glycerol with the overnight culture and store in 1 mL aliquots at -80 °C. This is the TG1 stock, which we use to prepare log-phase TG1 culture when needed.
- d. To 25 mL of 2× YT in a 125 mL disposable baffled Erlenmeyer flask, add enough TG1 stock (step D1c) so that so that the starting OD<sub>600</sub> is approximately 0.05. See **General note 12**.
- e. Incubate the TG1 culture at 37 °C shaking at 200 rpm until the OD<sub>600</sub> reaches 0.4–0.6. This takes approximately 1.5–2 h. This is the log-phase TG1, which should be kept on ice and used within 2 h.

2. Prepare the CM13K helper phage

**Caution:** Phage can infect all bacteria. If working in a protein expression lab, care must be taken not to release phage into communal areas. Wipe down surfaces with 5% Chemgene followed by 70% EtOH if the surface is susceptible to corrosion. The use of filtered tips is encouraged to reduce aerosols. All contaminated consumables must be decontaminated with 2% Virkon and autoclaved. If possible, use disposable Erlenmeyer flasks for all cultures outlined in steps D1–D5.

*Note: Work in a safety cabinet.*

- a. Inoculate 20 mL of 2× YT with 200 µL of log-phase TG1 (step D1e) in a 125 mL disposable baffled Erlenmeyer flask. Add 1 µL of previously prepared  $1 \times 10^{13}$  PFU/mL CM13K helper phage glycerol stock (step D2h) and incubate at 37 °C for 4 h shaking at 250 rpm.
- b. In two 1 L disposable baffled Erlenmeyer flasks, add 250 mL of 2× YT containing 25 µg/mL kanamycin and 10 mL of the prepared culture. Incubate overnight at 25 °C shaking at 250 rpm.
- c. Pellet the cells at 4,500× g for 10 min at 4 °C and transfer approximately 170 mL of the phage containing supernatant into three 250 mL PPCO centrifuge tubes. Add approximately 35 mL of ice-cold PEG/NaCl to each tube, invert, and swirl gently to mix the PEG/NaCl solution with the phage supernatant. Incubate on ice for 1 h.
- d. Centrifuge the tubes at 12,000× g for 15 min at 4 °C and resuspend each of the three pellets in 4 mL of ice-cold sterile PBS. Combine the resuspended pellets into a single 25 mL high-speed PPCO centrifuge tube and add 2.5 mL of ice-cold PEG/NaCl. Swirl and invert gently to mix the PEG/NaCl solution with the phage supernatant. Incubate on ice for 30 min.
- e. Centrifuge at 12,000× g for 10 min at 4 °C and resuspend the pellet in 12 mL of ice-cold sterile PBS. Keep at 4 °C until the phage concentration, expressed as plaque forming units per millilitre (PFU/mL), has been determined either spectrophotometrically or by plaque assay, both of which should yield similar values.
- f. For the spectrophotometric determination of phage concentration, in a 1.5 mL microcentrifuge tube, add 998 µL of PBS and 2 µL of phage (step D2e), transfer to a quartz cuvette, and measure the absorbance at 268 nm. Calculate the PFU/mL considering that an A<sub>268</sub> of 1.0 is equivalent to  $5 \times 10^{12}$  PFU/mL [11].
- g. For the determination of phage concentration using the plaque assay, prepare tenfold serial dilutions

( $10^{-1}$  to  $10^{-25}$ ) of CM13K helper phage (step D2e) in log-phase TG1 (step D1e) and incubate at RT for 5 min. Plate 100  $\mu$ L of  $10^{-15}$  to  $10^{-25}$  dilutions onto 1% agar LB plates without antibiotics. Add 3 mL of 0.7% agar LB (at 55 °C) per plate and swirl over the entire plate surface. Allow the top agar to solidify and incubate overnight at 37 °C. Calculate the PFU/mL using the following formula: number of plaques (clear areas)  $\times$  10 (to get to mL)  $\times$  dilution.

- h. Dilute the CM13K helper phage (step D2e) to  $2 \times 10^{13}$  PFU/mL using ice-cold sterile PBS, add an equal volume of ice-cold 50% glycerol, and mix gently to achieve a final concentration of  $1 \times 10^{13}$  PFU/mL. Store in 100  $\mu$ L aliquots in 500  $\mu$ L microcentrifuge tubes. This is the CM13K helper phage stock. Each time a new 100  $\mu$ L aliquot is thawed, sub-aliquot 12  $\mu$ L into 9  $\times$  PCR tubes to avoid freeze thawing of the original 100  $\mu$ L aliquot.

### 3. Panning first round

*Note: Work at a lab bench.*

- a. In a 2 mL microcentrifuge tube, add 500  $\mu$ L of StartingBlock buffer and 500  $\mu$ L of isolated library phage (step D3j). Wrap parafilm around the microcentrifuge tube lid to avoid any accidental spillage and incubate for 30 min at RT using a sample mixer.
- b. Add an appropriate volume of biotinylated antigen to achieve 50 nM final concentration in 1.5 mL (which is achieved by step D4d). See **General note 1**. Wrap parafilm around the microcentrifuge tube lid to avoid any accidental spillage and incubate for 1 h at RT using a sample mixer.
- c. At this point, one can start growing the log-phase TG1 cells (steps D1d–e), which would ensure that the cells are at the required OD<sub>600</sub> once the phage has been eluted (step D4h).
- d. Thirty minutes into the incubation of phage with antigen (step D4b), in a new 2 mL microcentrifuge tube, add 100  $\mu$ L of streptavidin Dynabeads™ M-280. Wash the beads twice with 500  $\mu$ L of PBS and then add 500  $\mu$ L of StartingBlock buffer. Incubate for 30 min at RT using a sample mixer.
- e. Transfer the 1 mL of blocked phage with antigen (step D4b) to the blocked beads (step D4d). The total volume is 1.5 mL, and the final antigen concentration is now 50 nM. Wrap parafilm around the microcentrifuge tube lid to avoid any accidental spillage and incubate for 15 min at RT using a sample mixer.
- f. Place the microcentrifuge tube on a DynaMag™-2 magnet. Once the beads have been pulled to the side, discard the unbound phage in the supernatant. Wash away loosely bound phage by resuspending the beads in 500  $\mu$ L of PBST and then using the magnet to pull the beads to the side. Repeat this washing process a total of six times with PBST and then once with PBS.
- g. Add 500  $\mu$ L of 250  $\mu$ g/mL trypsin to the beads. Wrap parafilm around the microcentrifuge tube lid to avoid any accidental spillage and incubate for 30 min at RT using a sample mixer.
- h. Collect the eluted phage from the first pan in the supernatant after using the DynaMag™-2 magnet to pellet the beads.

*Note: Work in a safety cabinet.*

- i. In a 50 mL conical centrifuge tube, add 500  $\mu$ L of eluted phage from the first pan (step D4h) to 10 mL of log-phase TG1 cells (step D4c) and incubate at 37 °C for 30 min without agitation.
- j. Pellet the cells at 2,800 $\times$  g for 10 min at 18 °C and resuspend in 1 mL of 2 $\times$  YT. Spread the entire volume until dry on one bioassay dish containing 1% agar LB with 100  $\mu$ g/mL ampicillin and incubate overnight at 37 °C.
- k. Resuspend the lawn of colonies in 6 mL of 2 $\times$  YT (containing a final concentration of 25% glycerol) and store in 1 mL aliquots in cryotubes at -80 °C. This is the pan 1 sub-library stock.

**Table 11. Experimental timetable for the identification of nanobody binders through panning, ELISA, and small-scale expression**

	Monday	Tuesday	Wednesday	Thursday	Friday
Week 1	Amplify VHH library (step D3)	Recover library phage (step D3) First panning (step D4)	Recover lawn of colonies (step D4) Amplify sub-library (step D5)	Recovery of sub-library phage (step D5) Second panning (step D6)	Take plates out of the incubator and place in fridge.
Week 2	Pick 93 colonies and grow overnight (step D7)	Prepare master plates (step D7) Start culture for anti-M13 ELISA (step D8)	Anti-M13 ELISA (step D8)	PCR, clean up, sequencing (step D9)	
Week 3	Overnight culture for small scale (step E1)	Small-scale culture to OD and induce with IPTG (step E1)	Periplasm isolation (step E1) Ni-NTA purification (step E2) Desalting (step E3)	Titration ELISA (step E4) Off-rate determination (step E5)	

4. Amplify and recover sub-library

*Note: Work in a safety cabinet.*

- a. To 50 mL of 2× YT containing 100 µg/mL ampicillin in a 250 mL disposable baffled Erlenmeyer flask, add enough pan 1 sub-library stock (step D4k) so that the starting OD<sub>600</sub> is approximately 0.05. See **General note 13**.
- b. Incubate the culture at 37 °C shaking at 200 rpm until the OD<sub>600</sub> reaches 0.4–0.6. This takes approximately 2 h. See **General note 14**.
- c. Transfer 10 mL of this culture into a 50 mL conical centrifuge tube, add 10 µL of the CM13K helper phage stock (step D2h), and mix by gentle swirling. Incubate at 37 °C for 1 h without agitation.
- d. Pellet the cells at 2,800× g for 10 min at 18 °C. Resuspend the pellet in 50 mL of 2× YT containing 100 µg/mL ampicillin and 25 µg/mL kanamycin and transfer to a 250 mL disposable baffled Erlenmeyer flask. Incubate overnight at 25 °C shaking at 250 rpm.
- e. Transfer the overnight culture into a 50 mL conical centrifuge tube and pellet the cells at 3,200× g for 10 min at 4 °C.
- f. Pour 40 mL of the phage containing supernatant into a new 50 mL conical centrifuge tube and add 10 mL of ice-cold PEG/NaCl. Invert several times and incubate on ice for 1 h.
- g. Pellet the precipitated phage by centrifugation at 3,200× g for 10 min at 4 °C. Resuspend the phage in 1 mL of ice-cold sterile PBS and transfer to a 2 mL microcentrifuge tube.
- h. After centrifugation at 20,000× g for 1 min at 4 °C, transfer the supernatant containing phage into a new 2 mL microcentrifuge tube. Add 250 µL of ice-cold PEG/NaCl and invert the microcentrifuge tube until a homogeneous white suspension appears. Incubate on ice for 30 min.
- i. Pellet the precipitated phage by centrifugation at 20,000× g for 15 min at 4 °C. Resuspend the pelleted phage in 1 mL of ice-cold sterile PBS and transfer to a 2 mL microcentrifuge tube.
- j. After a final centrifugation at 20,000× g for 1 min at 4 °C, transfer the supernatant containing phage into a new 2 mL microcentrifuge tube. This is the isolated pan 1 phage; it should be kept either on ice or at 4 °C and is stable for a month.

5. Panning second round

*Note: Work in a lab bench.*

- a. In a 2 mL microcentrifuge tube, add 500 µL of StartingBlock buffer and 500 µL of isolated pan 1 phage (step D5j). Wrap parafilm around the microcentrifuge tube lid to avoid any accidental spillage and incubate for 30 min at RT using a sample mixer.
- b. Add an appropriate volume of biotinylated antigen to achieve 5 nM final concentration in 1.5 mL (which is achieved by step D6d). See **General note 1**. Wrap parafilm around the microcentrifuge tube lid to avoid any accidental spillage and incubate for 1 h at RT using a sample mixer.
- c. At this point, one can start growing the log-phase TG1 cells (steps D1d–e), which would ensure that the cells are at the required OD<sub>600</sub> once the phage has been eluted (step D6h).
- d. Thirty minutes into the incubation of phage with antigen (step D6b), in a new 2 mL microcentrifuge tube, add 10 µL of streptavidin Dynabeads™ M-280. Wash the beads twice with 500 µL of PBS and then add 500 µL of StartingBlock buffer. Incubate for 30 min at RT using a sample mixer.
- e. Transfer the 1 mL of blocked phage with antigen (step D6b) to the blocked beads (step D6d). The total volume is now 1.5 mL, and the final antigen concentration is now 5 nM. Wrap parafilm around the microcentrifuge tube lid to avoid any accidental spillage and incubate for 15 min at RT using a sample mixer.
- f. Place the microcentrifuge tube on a DynaMag™-2 magnet. Once the beads have been pulled to the side, discard the unbound phage in the supernatant. Wash away loosely bound phage by resuspending the beads in 500 µL of PBST and then using the magnet to pull the beads to the side. Repeat this washing process a total of six times with PBST and then once with PBS.
- g. Add 500 µL of 250 µg/mL trypsin to the beads. Wrap parafilm around the microcentrifuge tube lid to avoid any accidental spillage and incubate for 30 min at RT using a sample mixer.
- h. Collect the eluted phage from the second pan in the supernatant after using the DynaMag™-2 magnet to pellet the beads.

*Note: Work in a safety cabinet.*

- i. To obtain individual colonies for the master plate (step D7), titrate the eluted second pan phage (step D6h) using tenfold serial dilutions from  $10^{-1}$  to  $10^{-12}$  with log-phase TG1 cells (step D6c). After incubation at 37 °C for 15 min, spread 50  $\mu$ L of the  $10^{-4}$  to  $10^{-12}$  dilutions on 1% agar LB plates containing 100  $\mu$ g/mL ampicillin until dry and incubate overnight at 37 °C. Include a TG1 only control. See **General note 15 and Troubleshooting 1**.
  - j. Select the dilution that has approximately 100–300 individual colonies.
6. Preparation of master plates

*Note: Work in a safety cabinet.*

- a. In a 1 mL deep-well block, add 180  $\mu$ L/well of 2 $\times$  YT containing 100  $\mu$ g/mL ampicillin. Pick 93 individual clones from selected plate of titrated pan 2 phage in TG1 (step D6j) with 93  $\times$  10  $\mu$ L tips and place into each filled well. Remove the tips using a multichannel pipette. Leave wells F12, G12, and H12 without a colony, which serve as our controls.
  - b. Cover the plate with cell culture adhesive seal and incubate overnight at 37 °C shaking at 600 rpm. This overnight culture will be used to prepare two stock plates, the PCR template (step D7c) and the master plate (step D7d). Prepare the PCR plate first followed by the master plate.
  - c. Transfer 20  $\mu$ L/well of overnight culture to a skirted PCR plate. Cover the plate with PCR foil seal and store at -80 °C. This is the PCR template plate.
  - d. To the remaining overnight culture, add 150  $\mu$ L/well of 50% glycerol and mix by pipetting. Cover the plate with PCR foil seal and store at -80 °C. This is the master plate.
7. Anti-M13 phage ELISA

*Note: Work in a safety cabinet.*

- a. In a 1 mL deep-well plate already filled with 150  $\mu$ L/well of 2 $\times$  YT, add 10  $\mu$ L/well of the master plate (step D7d). Cover with cell culture adhesive seal and incubate for 3 h at 37 °C shaking at 600 rpm.
- b. Prepare the CM13K helper phage mixture by adding 2  $\mu$ L of CM13K helper phage stock (step D2h) to 10 mL of 2 $\times$  YT. Add 30  $\mu$ L/well of the prepared CM13K helper phage mixture and incubate for 1 h at 37 °C shaking at 600 rpm.
- c. Pellet the cells at 2,800 $\times$  g for 10 min at 18 °C and resuspend each well in 400  $\mu$ L of 2 $\times$  YT containing 100  $\mu$ g/mL ampicillin and 25  $\mu$ g/mL kanamycin.
- d. Cover the plate with cell culture adhesive seal and incubate overnight at 25 °C shaking at 600 rpm.
- e. After centrifugation at 4,000 $\times$  g for 15 min at 4 °C, transfer the phage containing supernatant to a new 0.5 mL deep-well plate and cover with adhesive film. Store at 4 °C; this is stable for a month. This is the culture that is used for the anti-M13 ELISA (step D8l).

*Note: Work in a lab bench.*

- f. Coat a 96-well ELISA microplate with 100  $\mu$ L/well of 10  $\mu$ g/mL neutrAvidin™ diluted in PBS and incubate overnight at 4 °C.
- g. Using a plate washer, wash the plate five times with 300  $\mu$ L/well of PBST.
- h. Add 100  $\mu$ L/well of 50 nM biotinylated target protein diluted in PBS and incubate for 1 h at RT on a microplate shaker. See **General note 1**.
- i. Using a plate washer, wash the plate five times with 300  $\mu$ L/well of PBST.
- j. Add 250  $\mu$ L/well of blocking solution and incubate for 1 h at RT on a microplate shaker.
- k. Using a plate washer, wash the plate five times with 300  $\mu$ L/well of PBST.
- l. Add 50  $\mu$ L/well of blocking buffer followed by 50  $\mu$ L/well of the anti-M13 ELISA culture (step D8e). In the control wells F12, G12, and H12, add 50  $\mu$ L/well of blocking buffer and 50  $\mu$ L/well of PBS. Incubate for 1 h at RT on a microplate shaker.
- m. Using a plate washer, wash the plate five times with 300  $\mu$ L/well of PBST.
- n. Add 100  $\mu$ L/well of anti-M13-HRP diluted in 0.1% BSA-PBS (at a final concentration of 0.25  $\mu$ g/mL) and incubate for 1 h at RT on a microplate shaker.
- o. Using a plate washer, wash the plate five times with 300  $\mu$ L/well of PBST.
- p. Add 100  $\mu$ L/well of ABTS substrate, which is prepared by mixing solution A and solution B in a 1:1

- ratio. Protect from light and measure the absorbance at 405nm. See **General note 2**.
- q. Wells with absorbances above 1 are considered hits. See **Troubleshooting 2 and 3**.
8. PCR and sequencing of ELISA hits
- a. Prepare the following PCR reaction master mixture in a 5 mL tube and set up 96 PCR reactions in a 96-well rigid PCR plate (Table 12). See **General note 17**. Cover the plate with PCR foil seal.

**Table 12. PCR reaction composition**

Reagent	Final concentration	Volume for single reaction	Volume for master mixture
Nuclease-free water	n/a	18.25 $\mu$ L	1,825 $\mu$ L
10 $\mu$ M PhD seq Fwd primer	0.2 $\mu$ M	0.5 $\mu$ L	50 $\mu$ L
10 $\mu$ M PhD seq Rev primer	0.2 $\mu$ M	0.5 $\mu$ L	50 $\mu$ L
10 $\times$ Taq buffer	1 $\times$	2.5 $\mu$ L	250 $\mu$ L
25 mM dNTP	0.25 mM	0.25 $\mu$ L	50 $\mu$ L
Taq polymerase	1 unit	0.5 $\mu$ L	50 $\mu$ L
PCR template (step D7c)		2.5 $\mu$ L	-
Total	n/a	25 $\mu$ L	2275 $\mu$ L

- b. Set up the thermal cycler as below and perform the PCR (Table 13).

**Table 13. PCR cycling conditions**

Step	Temperature ( $^{\circ}$ C)	Duration	Number of cycles
Denaturation	95	7 min s	1
Annealing	95	15 s	
Extension	55	30 s	35
Final extension	68	1 min 40 s	
Hold	68	5 min	1
Denaturation	12	Infinite hold	-

- c. Add 45  $\mu$ L of Highprep PCR magnetic beads/well and mix by pipette. Place on DynaMag<sup>TM</sup> 96-side magnet for 1 min or until the solution is clear and the beads have been pulled against the wall of the PCR tube and pipette off the supernatant.
- d. Add 100  $\mu$ L of 70% EtOH and mix by pipette. Place on DynaMag<sup>TM</sup> 96-side magnet for 1 min or until the solution is clear and the beads have been pulled against the wall of the PCR tube and pipette off the supernatant. Repeat this washing step a total of three times.
- e. Leave beads to dry at RT for 10 min and then add 25  $\mu$ L of EB buffer (from the QIAprep Spin Miniprep Kit) to the beads and mix by pipette. Place on DynaMag<sup>TM</sup> 96-side magnet for 1 min or until the solution is clear and the beads have been pulled against the wall of the PCR tube and transfer 20  $\mu$ L of the purified PCR product containing supernatant to a new 96-well skirted PCR plate.
- f. If a nanophotometer that can measure multiple samples at a time (such as the Implen<sup>TM</sup> NanoPhotometer<sup>®</sup> N120) is available, measure the A<sub>280</sub> of each well. If a single sample nanophotometer is available, measure the A<sub>280</sub> of a random selection of clones. Expect around 15–25 ng/ $\mu$ L. Alternatively, 5  $\mu$ L of purified PCR product of a random selection of clones can be mixed with 5  $\mu$ L of 2 $\times$  DNA gel loading buffer and run alongside a lane of 1 $\times$  GeneRuler 1 kb DNA ladder on a 1% agarose gel containing 1 $\times$  SYBR<sup>TM</sup> Safe DNA gel stain. Electrophorese in 1 $\times$  TBE buffer at 80 V for 40 min. Image the gel and confirm the presence of amplified ~500 bp VHH genes.
- g. Cover the plate with PCR foil seal and send for sequencing along with the PhD seq Fwd primer. See **General note 18**.



## E. Small-scale expression and purification

The identification of nanobody clones that are expressed and purified well is achieved by performing a small-scale expression. The addition of a desalting step after purification allows for the confirmation of nanobody binding to the antigen using a titration ELISA or biolayer interferometry (BLI). The off-rate ( $k_{dis}$ ) of the nanobodies can be determined as well as preliminary epitope binning using BLI. These factors aid in the selection of which nanobodies will be taken forward for large scale expression.

### 1. Small-scale expression

*Note: Work in a safety cabinet.*

- In a 1 mL deep-well plate already filled with 180  $\mu$ L of 2 $\times$  YT containing 100  $\mu$ g/mL ampicillin per well, add 10  $\mu$ L of selected clones based on sequencing results from the master plate (step D7d). Cover with cell culture adhesive seal and incubate overnight at 37  $^{\circ}$ C shaking at 600 rpm.
- Add 10  $\mu$ L of the overnight culture of each clone to a 24-well plate already filled with 4 mL/well of terrific broth containing 100  $\mu$ g/mL ampicillin, 0.1% glucose, and 2 mM  $MgCl_2$ . Cover with cell culture adhesive seal and incubate for 4 h at 37  $^{\circ}$ C, shaking at 600 rpm.
- To each well, add 4  $\mu$ L of 1 M IPTG (1 mM final concentration), cover with cell culture adhesive seal, and incubate overnight at 25  $^{\circ}$ C shaking at 600 rpm.

*Note: Work in a lab bench.*

- Pellet the cells at 4,000 $\times$  g for 15 min at 4  $^{\circ}$ C and gently resuspend the pellet in 300  $\mu$ L of 1 mg/mL polymyxin B sulfate. Cover plate with adhesive film and agitate for 1 h at RT using a microplate shaker.
- Pellet the cells at 4,000 $\times$  g for 10 min at 4  $^{\circ}$ C and transfer the periplasm containing supernatant to a new 1.5 mL microcentrifuge tube.
- Mix 10  $\mu$ L of the isolated periplasm with 10  $\mu$ L of 2 $\times$  Laemmli sample buffer, incubate at 95  $^{\circ}$ C for 5 min, and run on a NuPage 4%–12% bis tris precast gel in 1 $\times$  MES buffer alongside Mark12<sup>TM</sup> unstained standard at 200 V for 40 min.
- Stain with InstantBlue<sup>®</sup> Coomassie protein stain. Expect prominent bands at approximately 14 kDa. Refer to the results from Expasy Protparam (Data analysis, step 1d) for the exact size of each clone.

### 2. Affinity purification

*Note: Work in a lab bench.*

- Purification of the his-tagged nanobodies using Ni-NTA spin columns is performed according to the manufacturer's instructions. Briefly, add 600  $\mu$ L of equilibration buffer to the Ni-NTA spin column and centrifuge at 2,900 rpm (800 $\times$  g) for 2 min at 4  $^{\circ}$ C.
- Add the periplasm sample to the Ni-NTA spin column and centrifuge at 1,600 rpm (200 $\times$  g) for 5 min at 4  $^{\circ}$ C.
- Wash the Ni-NTA spin column three times by the addition of 600  $\mu$ L of equilibration buffer and centrifuging at 2,900 rpm (800 $\times$  g) for 2 min at 4  $^{\circ}$ C.
- Transfer the Ni-NTA spin column into a new 2 mL microcentrifuge tube and add 100  $\mu$ L of elution buffer. Incubate at RT for 2 min before centrifugation at 2,900 rpm (800 $\times$  g) for 2 min at 4  $^{\circ}$ C. The flowthrough contains the purified nanobody.
- Mix 5  $\mu$ L of the purified nanobody with 5  $\mu$ L of 2 $\times$  Laemmli sample buffer, incubate at 95  $^{\circ}$ C for 5 min, and run on a NuPage 4%–12% bis tris precast gel in 1 $\times$  MES buffer alongside Mark12<sup>TM</sup> unstained standard at 200 V for 40 min.
- Stain with InstantBlue<sup>®</sup> Coomassie protein stain. Expect bands at approximately 14 kDa. Refer to the results from Expasy Protparam (Data analysis, step 1d) for the exact size of each clone.

### 3. Desalting

*Note: Work in a lab bench.*

- Desalting of the Ni-NTA purified nanobodies using the Zeba<sup>TM</sup> Spin Desalting Columns is performed according to the manufacturer's instructions. Briefly, break off the bottom cap of the desalting columns, loosen top cap, and place spin column in a new 2 mL microcentrifuge tube. Centrifuge at




- 1,500× g for 1 min at 4 °C.
  - b. Make a mark on the side of the desalting column where the resin is slanted upwards. Ensure that this mark is facing outwards in all centrifugation steps.
  - c. Wash the resin three times by gently adding 300 µL of PBS (or buffer of choice) to the resin bed and centrifuging at 1,500× g for 1 min at 4 °C.
  - d. Transfer the desalting column to a new 1.5 mL microcentrifuge tube and gently add 100 µL of purified nanobody (step E2d) to the centre of the resin bed. After the sample has been fully absorbed, gently add 15 µL of PBS (or buffer of choice) to the centre of the resin bed. Centrifuge at 1,500× g for 2 min at 4 °C and collect the desalted nanobody in the flowthrough.
  - e. Measure the  $A_{280}$  of the desalted nanobody and determine the molar concentration ( $\mu\text{M} = A_{280} / \epsilon \times 10^6$ ) as well as the concentration ( $\mu\text{g/mL} = \mu\text{M} \times \text{kDa}$ ) using the kDa and extinction coefficient ( $\epsilon$ ) obtained from ExPASy ProtParam for the individual nanobodies (Data analysis, step 1d).
4. Titration ELISA
- If an Octet R8 (or similar instrument) is not available, a titration ELISA will suffice to confirm binding of the purified nanobodies to the target protein.
- Note: Work in a lab bench.*
- a. Coat a suitable number of wells of a 96-well ELISA microplate with 100 µL/well of 10 µg/mL neutrAvidin™ diluted in PBS and incubate overnight at 4 °C.
  - b. Using a plate washer, wash the plate five times with 300 µL/well of PBST.
  - c. Add 100 µL/well of 50 nM biotinylated target protein diluted in PBS and incubate for 1 h at RT on a microplate shaker. See **General note 1**.
  - d. Using a plate washer, wash the plate five times with 300 µL/well of PBST.
  - e. Add 250 µL/well of blocking solution and incubate for 1 h at RT on a microplate shaker.
  - f. Using a plate washer, wash the plate five times with 300 µL/well of PBST.
  - g. Add 200 µL of 25 µg/mL desalted nanobody (step E3d) diluted in blocking solution to the first well. Prepare a serial dilution in the plate, i.e., add 100 µL of blocking solution to the second well and add 100 µL of the solution from the first well. Prepare a total of seven dilutions (25–0.39 µg/mL), leaving the last well with blocking solution only. Incubate for 1 h at RT on a microplate shaker.
  - h. Using a plate washer, wash the plate five times with 300 µL/well of PBST.
  - i. Add 100 µL/well of anti-camelid VHH-HRP diluted 1:5,000 in 0.1% BSA-PBS and incubate for 1 h at RT on a microplate shaker.
  - j. Using a plate washer, wash the plate five times with 300 µL/well of PBST.
  - k. Add 100 µL/well of ABTS substrate, which is prepared by mixing solution A and solution B in a 1:1 ratio. Protect from light and measure the absorbance at 405 nm. See **General note 2**.
5. Off-rate ( $k_{\text{dis}}$ ) determination using Biolayer interferometry (BLI)
- Note: Work in a lab bench.*
- a. Preheat the Octet R8 at 25 °C overnight or at least 3 h before starting the experiment.
  - b. Typical starting concentrations of antigen and nanobody are 50 nM and 100 nM, respectively, prepared in octet dilution buffer. See **General note 19**.
  - c. Prepare all the dilutions in a 96-well black polypropylene plate. See **General note 20**.
  - d. A single replicate per nanobody is sufficient at this stage. Include a no nanobody control to rule out any non-specific interactions with the loaded antigen and to serve as the background.
  - e. The experimental run conditions using SA biosensors are as follows:  
 1,000 rpm shaking for all incubation steps.  
 Plate definition: Octet dilution buffer to be used in all baseline steps/wells.  
 Assay definition: Baseline 600 s into octet dilution buffer.  
     : Loading 300 s into antigen diluted in octet dilution buffer.  
     : Baseline 30 s in octet dilution buffer.  
     : Baseline 60 s in octet dilution buffer.  
     : Association 300 s in desalted nanobody diluted in octet dilution buffer (watch this step, time might need to be extended as we want the curves to reach a plateau).  
     : Dissociation 300 s into column where baseline 60 s was measured (this step should

be the same duration as the association step).

F. Subcloning of the nanobody into the bacterial and mammalian vectors

We have developed two vector suites for the functionalisation of nanobodies in both bacterial and mammalian expression systems. They have been designed in such a way that any nanobody that has been identified from a library that has been constructed as detailed in section C can be cloned into them using a common forward primer and tag-specific reverse primers as well as using the same restriction sites (Figure 8). The characteristics, tags for functionalisation, and primers for each vector in the bacterial and mammalian toolboxes are detailed in Table 14.

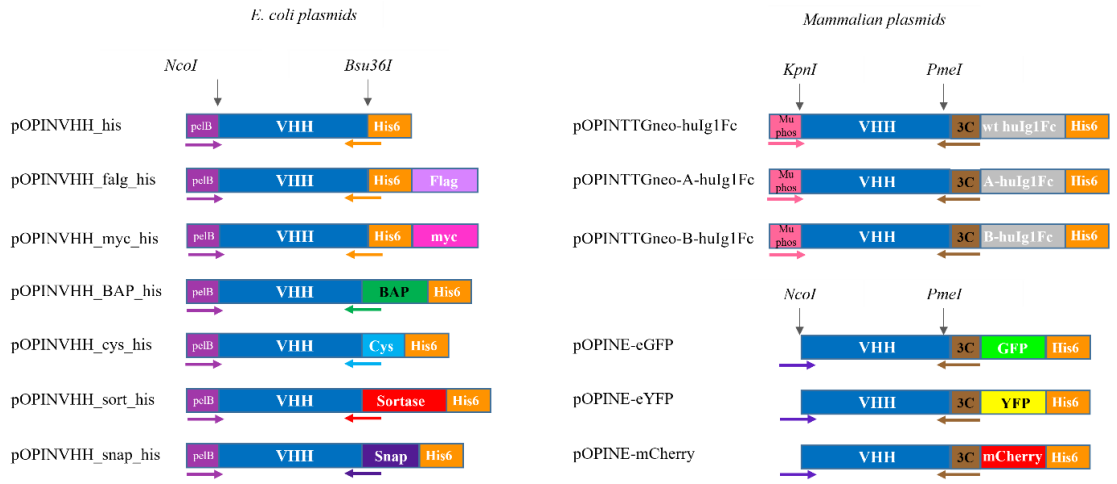

**Figure 8. Details of the primers, restriction enzymes, and tags present in each of the vectors in the bacterial and mammalian toolboxes.** Abbreviations used: VHH, nanobody; BAP, biotin acceptor peptide/Avi-Tag™; wt, wild-type; huIg1Fc, human Ig1 Fc.

Depending on the desired application, the selected nanobody is subcloned into a suitable expression vector. Below we detail the subcloning of a hit nanobody into the pOPINVHH\_His vector as an example. The same process is followed to subclone the same nanobody into a mammalian vector for mammalian cell expression.

1. Restriction digestion of the pOPINVHH\_His vector

*Note: Work at a lab bench.*

- a. Prepare the following digestion mixture in a PCR tube (1 × 100 µL) (Table 15).

**Table 15. Restriction digestion reaction composition**

Reagent	Final concentration	Volume
Nuclease-free water	n/a	80.5 µL
10× rCutSmart™ buffer	1×	10 µL
NcoI	2.5 units	2.5 µL
Bsu36I	2.5 units	2.5 µL
400 ng/µL vector DNA	18 µg	4.5 µL
Total	n/a	100 µL

- b. Briefly vortex and centrifuge the reaction tube.
- c. Incubate at 37 °C for 1 h and then at 80 °C for 20 min.
- d. Add 4× volume of B3 buffer (400 µL) and purify using 1× Purelink column from the Purelink PCR purification kit according to the manufacturer's instructions. Elute using 50 µL of E1 buffer. Measure the A<sub>280</sub> and expect a concentration between 60 and 100 ng/µL.

**Table 14. Features of the vector in the *E. coli* and mammalian toolboxes.** The bp within the primer that is complementary to the vector is indicated in italics.

Expression system	Vector	Tag	Application	Restriction enzymes	Forward primer (5'–3')	Reverse primer (5'–3', reverse complement)
<i>E. coli</i>	pOPINVHH_his	His	Crystallography	NcoI and Bsu36I	<i>GCGGCCAGCCGGCC</i>	<i>GTGATGGTGGCCTGAG</i>
	pOPINVHH_flag_his	His and FLAG	Epitope detection		<i>ATGGCCAGGTGCAG</i>	<i>GAGACGGTGACCTGGG</i>
	pOPINVHH_myc_his	His and myc	Epitope detection		<i>CTGGTGGAG</i>	<i>TC</i>
	pOPINVHH_BAP_his	Avi-tag <sup>TM</sup> and His	Biotinylation			<i>ATCATTCAAGCCTGAGG</i>
						<i>AGACGGTGACCTGGGT</i>
						<i>C</i>
	pOPINVHH_cys_his	Cys and His	Chemical conjugation			<i>ATGGTGACAGCCTGAGG</i>
Mammalian						<i>AGACGGTGACCTGGGT</i>
						<i>C</i>
	pOPINVHH_sort_his	Sortase	Enzymatic labelling			<i>CGGCAGGCCGCCTGAG</i>
						<i>GAGACGGTGACCTGGG</i>
						<i>TC</i>
	pOPINVHH_snap_his	SNAP-tag <sup>®</sup>	Enzymatic labelling			<i>GTCCTTGTCGCCTGAGG</i>
						<i>AGACGGTGACCTGGGT</i>
Mammalian						<i>C</i>
	pOPINTTgneo-huIg1Fc	huIg1Fc and His	Fc fusions	KpnI and PmeI	<i>GCGTAGCTGAAACCGG</i>	<i>CAGAACTTCCAGTTTAG</i>
	pOPINTTgneo-huIg1FcA (Y407T)	huIg1Fc and His	Bispecific nanobody generation		<i>CCAGGTGCAGCTGGT</i>	<i>GGGAGACGGTGACCTG</i>
	pOPINTTgneo-huIg1FcB (T334Y)	huIg1Fc and His	Bispecific nanobody generation		<i>GGAG</i>	<i>GGTC</i>
	pOPINE-eGFP	GFP and His	Fluorescent fusions	NcoI and PmeI	<i>AGGAGATATACCATGC</i>	
	pOPINE-eYFP	YFP and His	Fluorescent fusions		<i>AGGTGCAGCTGGTGG</i>	
	pOPINE-mCherry	mCherry and His	Fluorescent fusions		<i>AG</i>	

2. Preparation of nanobody insert

*Note: Work in a safety cabinet.*

- a. From the master plate (step D7d), add 10  $\mu$ L of the selected nanobody clones to a 24-well plate already filled with 4 mL/well of 2 $\times$  YT containing 100  $\mu$ g/mL ampicillin, cover with cell culture adhesive seal, and incubate overnight at 37  $^{\circ}$ C, shaking at 600 rpm.

*Note: Work at a lab bench.*

- b. Use 2 mL of the overnight culture to isolate the plasmid using the QIAprep spin miniprep kit according to the manufacturer's protocol. Elute in a final volume of 50  $\mu$ L of EB. Measure the  $A_{280}$  and expect a concentration between 100 and 500 ng/ $\mu$ L.
- c. Send an aliquot for sequencing along with the PhD seq Fwd primer. If it matches the sequence obtained by PCR (step D9g), continue to the next step. See **General note 18**.
- d. Prepare the following reaction in a PCR tube (1  $\times$  25  $\mu$ L) (Table 16).

**Table 16. PCR reaction composition**

Reagent	Final concentration	Volume
Nuclease-free water	n/a	9 $\mu$ L
10 $\mu$ M Common Fwd primer	0.5 $\mu$ M	1.25 $\mu$ L
10 $\mu$ M His FLAG Rev primer	0.5 $\mu$ M	1.25 $\mu$ L
2 $\times$ Phusion flash PCR master mix	1 $\times$	12.5 $\mu$ L
10 ng/ $\mu$ L vector template (step F2b)	10 ng	1 $\mu$ L
Total	n/a	25 $\mu$ L

- e. Briefly vortex and centrifuge the reaction tube.
- f. Set up the thermal cycler as below and perform the PCR (Table 17).

**Table 17. PCR cycling conditions**

Step	Temperature ( $^{\circ}$ C)	Duration	Number of cycles
Denaturation	98	10 s	1
Annealing	98	1 s	
Extension	60	1 s	30
Final extension	72	5 s	
Hold	72	1 min	1
Denaturation	12	Infinite hold	-

- g. Add 5  $\mu$ L of 6 $\times$  DNA gel loading buffer to the reaction tube and run alongside a well of 1 $\times$  GeneRuler 1 kb DNA ladder on a 1% agarose gel containing 1 $\times$  SYBR<sup>TM</sup> safe DNA gel stain. Electrophorese in 1 $\times$  TBE buffer at 80 V for 40 min.
- h. Visualise the ~500 bp amplified nanobody band using a LED Blue transilluminator, excise from the gel using a scalpel knife, and gel extract using one spin column from the Nucleospin Gel and PCR clean-up kit according to the manufacturer's protocol. Elute using 50  $\mu$ L of NE, measure the  $A_{280}$ , and expect a concentration between 30 and 100 ng/ $\mu$ L.

**Caution:** Use Safe Imager<sup>TM</sup> viewing glasses when visualising the DNA during excision of the band from the agarose gel and use caution when using scalpel knives.

3. Ligation and transformation

*Note: Work at a lab bench.*

- a. Prepare the following reaction in a PCR tube (1  $\times$  10  $\mu$ L) (Table 18).

**Table 18. In-Fusion reaction composition**

Reagent	Final concentration	Volume
Nuclease-free water	n/a	5 $\mu$ L

5× CEII buffer	1×	2 µL
4 ng/µL insert (step F2h)	4 ng	1 µL
19.83 ng/µL vector DNA (step F1d)	19.83 ng	1 µL
Exnase II		1 µL
Total	n/a	10 µL

- b. Briefly vortex and centrifuge the reaction tube.
- c. Incubate at 42 °C for 30 min using a thermal cycler. Add 40 µL of TE buffer to terminate the reaction.  
*Note: Work in a safety cabinet.*
- d. Gently mix 5 µL of in-fusion reaction with 25 µL of Stellar™ competent cells in a sterile 2 mL microcentrifuge tube. Incubate on ice for 10 min.
- e. Heat shock at 42 °C for 45 s and incubate on ice for 2 min. Add 200 µL of LB medium and incubate for 1 h at 37 °C agitating at 600 rpm.
- f. Spread 50 µL of culture on a 1% agar LB plate containing 100 µg/mL ampicillin, 2 mM IPTG, and 40 µg/mL X-gal until dry. Incubate overnight at 37 °C. Expect many white colonies and few-to-no blue colonies.
- g. Transfer two white colonies from the transformation plate to two wells of a 24-well plate already filled with 4 mL/well of 2× YT containing 100 µg/mL ampicillin. Cover with cell culture adhesive seal and incubate overnight at 37 °C, shaking at 600 rpm.  
*Note: Work at a lab bench.*
- h. Use 2 mL of the overnight culture to isolate the plasmid using the QIAprep spin miniprep kit according to the manufacturer's protocol. Elute in a final volume of 50 µL of EB. Measure the A<sub>280</sub> and expect a concentration between 100 and 500 ng/µL.
- i. Send an aliquot for sequencing along with the PhD seq Fwd primer. If it matches the sequence obtained by PCR (step D9g), continue to the next step. See **General note 18**.
- j. Proceed to large-scale expression using WK6 *E. coli* or Expi293™ as detailed in Le Bas et al. [12] or use in imaging experiments using other mammalian cell lines (see section G) depending on the expression vector that was selected.

## G. Nanobody functionalisation and application

By using a selection of vectors from the *E. coli* and mammalian toolboxes, we demonstrate how nanobodies can be used as a reagent for pull downs using streptavidin and the pOPINVHH\_BAP vector, for imaging by chemical conjugation using the pOPINVHH\_Cys vector and by transfection with the pOPINE-GFP and pOPINE-YFP for in vivo and live cell imaging.


### 1. Confirmation of biotinylation of anti-GFP nanobody

The anti-GFP nanobody was cloned into the pOPINVHH\_BAP vector, expressed at large scale in BL21(DE3)-R3-pRARE2-BirA to allow for in vivo biotinylation, and purified using Ni-NTA and size exclusion chromatography as detailed in section F and in Le Bas et al. [12] with some modifications. Briefly, LB agar plates contained 34 µg/mL chloramphenicol, 50 µg/mL spectinomycin, and 50 µg/mL carbenicillin, cultures of terrific broth were supplemented with 50 µg/mL spectinomycin and 50 µg/mL carbenicillin, and a final concentration of 0.2 mM biotin was added alongside 0.1 mM IPTG. Below, we detail how we confirm that the anti-GFP nanobody has been biotinylated.

*Note: Work at a lab bench.*

- a. Add 2 µg of anti-GFP expressed in pOPINVHH\_BAP or pOPINVHH\_His and 1.5 µL of 50 µg/mL Streptavidin Alexa Fluor™ 488 conjugate to a 1.5 mL microcentrifuge tube and make up the volume to 10 µL with PBS. Incubate for 30 min at 37 °C.
- b. Add 10 µL of 2× Laemmli sample buffer, incubate at 95 °C for 5 min, and run on a NuPage 4%–12% bis tris precast gel in 1× MES buffer alongside Mark12™ unstained standard and BenchMark™ fluorescent protein standard at 200 V for 40 min.

- c. Image the in-gel fluorescence using the ProQ emerald 488 image acquisition preset on the ChemiDoc™ imaging system (blue Epi excitation, 532/28 emission filter). Stain with InstantBlue® Coomassie protein stain (Figure 9).

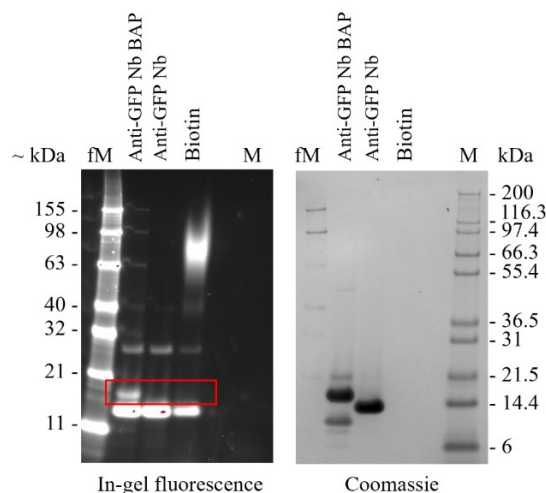

**Figure 9. Confirmation of nanobody biotinylation after expression in pOPINVHH\_BAP and co-expression with BirA by analysis on SDS-PAGE.** Samples were visualised by in-gel fluorescence and stained with InstantBlue® (Coomassie). Protein molecular weight markers in kDa were run in parallel.

## 2. Biotin-mediated pull downs

Below, we detail the utilisation of purified biotinylated anti-GFP nanobody in the purification of GFP from spiked *E. coli* lysate.

*Note: Work at a lab bench.*

- a. Wash 100  $\mu$ L of streptavidin Dynabeads™ M-280 twice with 1 mL of PBS.
- b. Add 100  $\mu$ L of 10  $\mu$ g biotinylated anti-GFP nanobody diluted in PBS to the prepared beads and incubate for 30 min at RT using a sample mixer.
- c. Collect the unbound sample and wash the beads three times with 1 mL of PBS.
- d. Add 10  $\mu$ g of pure GFP to *E. coli* lysate and make up to a final volume of 100  $\mu$ L using PBS. Add this to the washed beads and incubate for 30 min at RT using a sample mixer.
- e. Collect the unbound sample and wash the beads three times with 1 mL of PBS.
- f. Add 20  $\mu$ L of 2 $\times$  Laemmli sample buffer to the beads.
- g. To 10  $\mu$ L of each loaded and unbound sample, add 10  $\mu$ L of 2 $\times$  Laemmli sample buffer.
- h. Incubate all Laemmli-treated samples at 95 °C for 5 min and run on a NuPage 4%–12% bis tris precast gel in 1 $\times$  MES buffer alongside Mark12™ unstained standard at 200 V for 40 min.
- i. Stain with InstantBlue® Coomassie protein stain (Figure 12A).


## 3. Fluorophore conjugation

The anti-vimentin (VB3) nanobody was cloned into the pOPINVHH\_Cys vector, expressed at large scale in WK 6 *E. coli*, and purified using Ni-NTA and size exclusion chromatography as detailed in section F. Below, we detail how to prepare fluorescently labelled VB3 nanobody (VB3-AF647) and its application in *in vivo* imaging.

*Note: Work at a lab bench.*

- a. To 300  $\mu$ g of purified VB3 nanobody, add 1 mM TCEP pH 8.0 to make up to a final volume of 100  $\mu$ L and incubate for 20 min at 4 °C using a sample mixer.
- b. Add 1  $\mu$ L of Alexa Fluor™ 647 C<sub>2</sub>-maleimide to the VB3 nanobody solution and incubate for 1 h at 4 °C using a sample mixer.

- c. Remove any unbound Alexa Fluor™ 647 from the fluorescently labelled VB3 nanobody using a Zeba™ dye and biotin removal spin column as per the manufacturer's protocol. This sample is stored at 4 °C for two weeks protected from light but can be flash frozen and stored at -80 °C for an extended period.
- d. Mix 10 µL of the VB3-AF647 with 10 µL of 2× Laemmli sample buffer, incubate at 95 °C for 5 min, and run on a NuPage 4%–12% bis tris precast gel in 1× MES buffer alongside PageRuler™ prestained protein ladder at 200 V for 40 min.
- e. Image the in-gel fluorescence using the Oriole image acquisition preset on the ChemiDoc™ imaging system (UV excitation, 590/110 emission filter). Stain with InstantBlue® Coomassie protein stain (Figure 10).
- f. Measure the absorbance of the VB3-AF647 from 200 to 800 nm with baseline correction at 750 nm using the Nanodrop. Expect a peak at 280 and 651 nm, which correspond to VB3 nanobody and to Alexa Fluor™ 647, respectively (Figure 10).

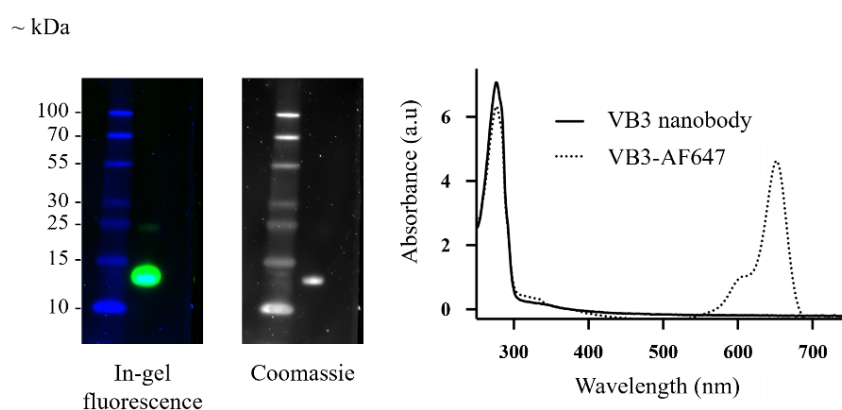

**Figure 10. Conjugation of Alexa Fluor™ 647 to VB3 nanobody (VB3-AF647) expressed in pOPINVHH\_Cys as determined by in-gel fluorescence (left), InstantBlue® staining (middle), and absorbance profile (right)**

#### 4. In vivo visualisation

*Note: Work in a safety cabinet.*

- a. Place a single coverslip into the wells of a 6-well culture plate. Seed 1 mL of 10,000 HeLa cells/well and culture in phenol red-free DMEM medium supplemented with 10% FBS, 1× GlutaMAX™, and 1× penicillin-streptomycin at 37 °C and 5% CO<sub>2</sub> in a CO<sub>2</sub> incubator.

*Note: Work at a lab bench.*

- b. After 24 h, wash the cells three times with 2 mL of PBS and fix with 2 mL of 4% paraformaldehyde for 20 min at RT.
- c. Wash the cells three times with 2 mL of PBS and block the cells with 2 mL of confocal blocking buffer for 1 h at RT.
- d. During the blocking step, prepare the fluorescently labelled VB3 nanobody solution for labelling vimentin in HeLa cells. In a 1.5 mL microcentrifuge tube, add 5 µL of VB3-AF647 (step G3c) to 1 mL of confocal dilution buffer PBS and mix gently by inversion.
- e. Remove the blocking buffer from the cells, add the prepared VB3-AF647 (step G4d), and incubate for 2 h at 4 °C using a gyratory rocker.
- f. Remove the VB3-AF647 solution and wash the cells three times with 2 mL of PBS.
- g. For each coverslip, add a small drop of Fluoroshield mounting media with DAPI onto a microscope slide. Using Dumont tweezers, gently remove the coverslip from the well and place cell side down onto the mounting media. Remove excess liquid by gently pressing the coverslip slide down towards



the microscope slide and secure the coverslip by adding nail varnish around the perimeter of the coverslip.

- h. The stained cells on the prepared slides are imaged using a Leica SP8 confocal microscope at 63× objective with type F immersion liquid. Use Ex405 and Em430–550nm for DAPI detection and Ex633 and Em650 for Alexa Fluor™ 647 detection.
5. Live cell imaging  
The VB3 nanobody was cloned into the pOPINE-eGFP and -eYFP vectors as detailed in section F. The produced vectors are referred to as pOPINE-VB3-eGFP and pOPINE-VB3-eYFP. Below, we detail how to transiently transfect HeLa cells with the pOPINE-VB3-eGFP and pOPINE-VB3-eYFP vectors for live-cell imaging of vimentin filaments.

*Note: Work in a safety cabinet.*

- a. In an 8-well chambered coverslip, seed 200 µL containing 25,000 HeLa cells/well and culture in phenol red-free DMEM medium supplemented with 10% FBS, 1× GlutaMAX™, and 1× penicillin-streptomycin at 37 °C and 5% CO<sub>2</sub> in a CO<sub>2</sub> incubator.
- b. After 24 h, prepare the transfection mixture as follows. In a 500 µL microcentrifuge tube, add 10 µL of serum-free DMEM medium and gently mix in 2 µL of 100 ng/µL pOPINE-VB3-eGFP and pOPINE-VB3-eYFP vectors. Add 0.5 µL of Fugene transfection reagent and incubate for 10 min at RT.
- c. Remove the old media from each of the wells in the 8-well microscope chambered coverslip and add 270 µL/well of DMEM medium supplemented with 10% FBS. Add each DNA-Fugene complex dropwise, by gently depressing the plunger of the pipette, into each well and incubate the chambered coverslip at 37 °C and 5% CO<sub>2</sub> in a CO<sub>2</sub> incubator for 48 h.
- d. Live-cell imaging is captured using a Leica SP8 confocal microscope at 63× objective with immersion liquid. Use Ex405 and Em430–550nm for DAPI detection and Ex488 and Em518-558 for GFP and YFP detection.


## Data analysis

1. Analysis of sequencing data (section D9)
  - a. Drag and drop the .seq files into [IMGT/V-QUEST](#). Make the following selections on the site: Species: Vicugna Pacos (alpaca); Receptor type or locus: Ig; choose to get the results in excel format. Submit. See **General note 21**.
  - b. There are multiple tabs in the produced Excel file. The fifth tab (AA-sequences) is the one that we will work in. Sort the sequences by productive, non-productive, and no results.
  - c. Then, sort the CDR3-IMGT of the productive sequences to visualise the CDR3 clusters. Highlighting the different CDR3 clusters in different colour aids in the selection of representative clones, which will be used for small-scale expression.
  - d. Open the .seq files of the selected clones in SnapGene to obtain the amino acid sequences. Copy the single letter amino acid sequence and paste into [Expasy - ProtParam tool](#). Delete the pelB sequence (MKYLLPTAAAGLLLLAAQPAMA), submit, and note the expected kDa, pI, and extinction coefficient (ε), assuming that all Cys residues are reduced.
2. Analysis of BLI data for off-rate determination (section E5)
  - a. In the *Preprocess data* tab: Subtract the background of the no nanobody well to all the other tested nanobody wells. Correct the Y axis to the baseline step and the interstep to the dissociation step. Apply the Savitzky-Golay filtering.
  - b. In the *Kinetic analysis* tab: Analyse both the association and dissociation steps. Use a 1:1 binding model. Use a local (individual) fitting. Fit the full association and dissociation steps and note the calculated off-rate ( $k_{dis}$ ).

## Validation of protocol

### Identification of nanobodies to the G protein of Nipah virus

The improved workflow for nanobody generation was exemplified by the identification of nanobodies to the receptor-binding G protein of the Nipah virus (NivG). The extracellular region of the NivG protein (residues 183–602) with N-terminal hexahistidine tag was expressed in Expi293<sup>TM</sup> cells and purified by a combination of immobilised metal affinity chromatography (IMAC) and size exclusion (Figure 11A). Comparison of pre- and post-immune sera in a seroconversion ELISA of biotin-tagged NivG confirmed that an immune response had been generated to the antigen (Figure 11B). Following two rounds of panning of the VHH library 13, 93 phage clones were picked and tested in a phage ELISA. The sequences of phage binders were amplified by PCR and sequenced and the translated VHHs were clustered by CDR3 sequence identity. A representative clone from each of the three major clusters was selected, expressed at small scale (4 mL cultures), and purified by IMAC (panel C). Binding to NivG was confirmed by titration ELISA (panel D). Analysis of binding off-rates by BLI confirmed that the three clones (A8, B7, and D5) bound to NivG (panel E) and were prioritised for large-scale expression.

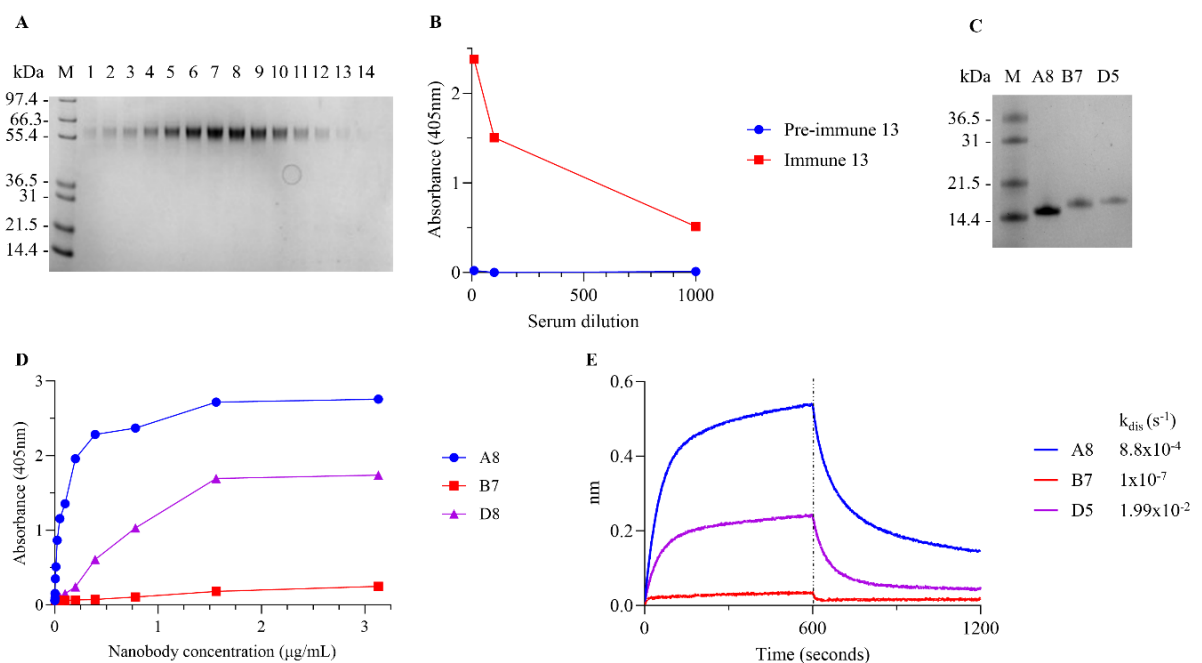


**Figure 11. Generation of NivG-binding nanobodies from an immunised llama.** (A) SDS-PAGE of fractions from gel filtration of purified NivG, (B) seroconversion ELISA of pre- and post-immune sera following immunisation with NivG, (C) SDS-PAGE of purified anti-NivG nanobodies from small-scale expression, (D) binding of purified nanobodies to NivG 183-602 by titration ELISA, (E) BLI traces of nanobody association to and dissociation from NivG 183-602; off-rates for the three nanobodies are shown.

### Functionalisation of nanobodies for purification and imaging

The functionalisation of nanobodies using three of the vectors presented in Table 1 and their application was exemplified using an anti-GFP [9] and an anti-vimentin (VB3) [10] nanobody. The anti-GFP nanobody was expressed in the pOPINVHH\_BAP vector, biotinylated *in vivo*, and used for pull-down of GFP from spiked *E. coli* lysate with magnetic-coupled streptavidin beads (Figure 12A). The anti-VB3 nanobody was expressed in *E. coli* using the pOPINVHH-cys vector, purified, and coupled to AlexaFluor<sup>TM</sup> 647 via a maleimide linker. HeLa cells were fixed and permeabilised and the vimentin filaments were stained with the anti-VB3 nanobody conjugate (Figure 12B). The same nanobody was genetically fused to either YFP or GFP by cloning into pOPINE-YFP and pOPINE-GFP vectors, respectively. Live HeLa cells were transfected with these vectors and imaged by confocal

microscopy to visualise vimentin filaments (Figure 12C).

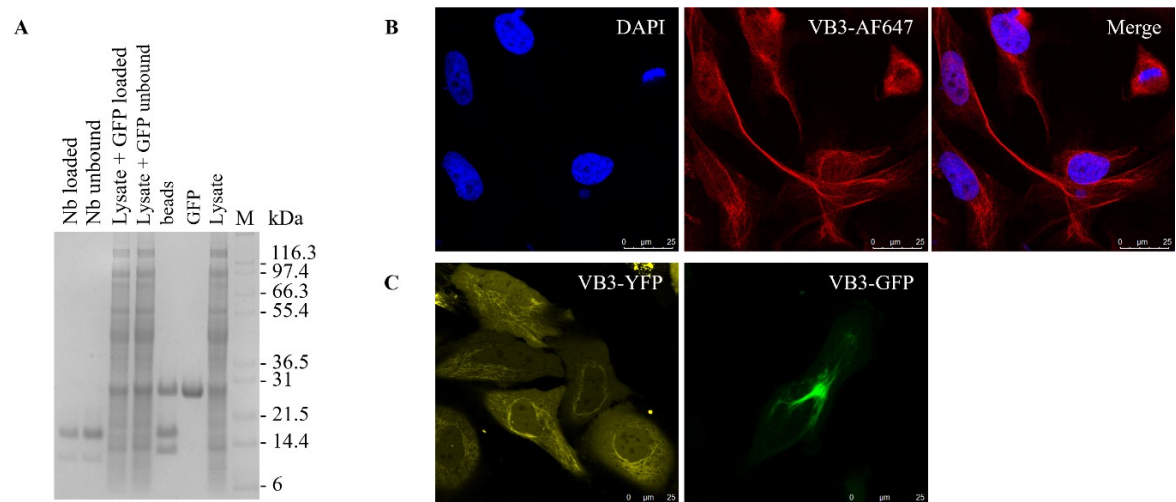

**Figure 12. Functionalisation of nanobodies.** (A) Pull-down of GFP spiked into an *E. coli* lysate by biotinylated anti-GFP nanobody coupled to streptavidin magnetic beads. Detection of vimentin in HeLa cells using (B) fluorescently labelled anti-vimentin nanobody (VB3-AF647) and (C) transiently transfected with VB3-eGFP and VB3-eYFP.

## General notes and troubleshooting

### General notes

1. Preparation of biotinylated proteins can be achieved in two ways: chemical conjugation or by in vitro biotinylation of a protein with the Avi-tag® sequence and co-expression with BirA.
2. Expect colour development within 5–15 min of substrate addition. If signals are low, the longest incubation before absorbance measurement is 1 h.
3. If we receive 170 mL of llama blood in 17 × 10 mL blood tubes without anticoagulant, the addition of an equal volume of PBS will result in a total of 340 mL of diluted blood. Fifteen millilitres of diluted blood are required for each 50 mL tube with 15 mL Histopaque®-1007 already added; therefore, a total of 22.6 tubes are required. As such, the 23rd tube will have less than 15 mL of diluted blood.
4. Holding the 50 mL conical centrifuge tube at a 45° angle allows for the gentle addition of the diluted blood over the Histopaque®-1007 to prevent the disturbance of the Histopaque®-1007 layer. During the process of gently layering the blood on top of the Histopaque®-1007 over the many 50 mL conical centrifuge tubes, the lower layer will turn red as the heavier red blood cells fall to the bottom of the tube. This is to be expected.
5. Mixing the trypan blue with the PMNCs too vigorously may damage the PMNCs, resulting in low viability values that may not be a true representation of the actual viability.
6. If no PCR workstation or dedicated RNase free space is available, keeping boxes of tips and bags microcentrifuge tubes that are used only for RNA work will be sufficient if gloves are changed regularly and RNaseZap™ is used.
7. RNA should be stored at -20 °C under ethanol where it is stable for extended periods of time. Once it is resuspended in water, the stability is reduced and, as such, for the best possible results, we recommend that cDNA synthesis is carried as soon as RNA has been resuspended in water (step C3i).
8. For the electroporations, expect ms values greater than 4.5.
9. Note that the estimated library sizes for the small-scale library will vary slightly from that of the scaled library.

10. An example on how to calculate how many microlitres of in-fusion reactions would be needed to make a library with a size of  $1 \times 10^6$ : We purify the 10  $\mu$ L in-fusion reaction using a Purelink column, eluting in 20  $\mu$ L nuclease-free water. We use 5  $\mu$ L of this for the electroporation, which translates to 2.5  $\mu$ L of the original in-fusion reaction being added to 30  $\mu$ L of competent TG1. Hence, if 2.5  $\mu$ L of in-fusion library yielded 373 CFU at 1:100 dilution, the small-scale library size would be  $1.492 \times 10^5$  CFU/mL. To create a  $1 \times 10^6$  CFU/mL library,  $6.7 \times 2.5$   $\mu$ L of in-fusion reactions would be required, which is just over 1  $\times 10$   $\mu$ L of in-fusion reaction. If enough VHH2 and digested pADL23c backbone has been prepared, we routinely use more in-fusion reactions than necessary, i.e., 10  $\times 10$   $\mu$ L, to ensure that we have the best chance of producing a  $1 \times 10^6$  CFU/mL library.
11. If there is a large discrepancy between the library size as calculated by titre and by plate density, use the plate density to guide your decision if the library is big enough. See Figure 7.
12. Begin by using 200  $\mu$ L of the TG1 stock when preparing the starting culture.
13. Begin by using 20  $\mu$ L of the isolated library phage and pan 1 sub-library stock when preparing the starting culture.
14. In some cases, it might take longer than 2 h for the culture to reach an OD<sub>600</sub> of 0.5. Continue culturing until the required OD<sub>600</sub> is reached. As long as the TG1 control plate in either step D3k or D6i are negative, you can be assured that this might just be variance in the library and/or sub-library.
15. A TG1-only control functions to show that the TG1s used have not been infected with phage prior to this titration step. See **Troubleshooting 1 and 2**.
16. Steps C8h–i show how diverse the library is. The titration of the amplified and recovered library (step D3l) gives an indication of how many phage particles we can make using our diverse library.
17. Adjust the volume of the prepared PCR master mixture according to the number of actual ELISA positive clones.
18. Volumes and concentrations of the DNA and primers might be supplier specific. Check the requirements of your local sequencing company.
19. Use enough antigen to load the biosensors to result in a 1 nm signal increase.
20. Polypropylene plates must be used as they do not bind protein.
21. The IMGT site can only process 48 sequences at a time. If more than 48 clones have been sequenced, combine the two generated IMGT Excel files together before sorting by V-DOMAIN functionality and CDR3-IMGT.

## Troubleshooting

Problem 1: Growth on TG1 control plates.

Possible cause: Indication of phage contamination of the *E. coli* before our isolated/eluted phage was added.

Solution: Prepare fresh log-phase TG1 from fresh aliquot of TG1 in the -80 °C. Repeat the titration of the phage in question.

Problem 2: Growth on TG1 control plates.

Possible cause: Ampicillin stock is old.

Solution: Using a new ampicillin aliquot, prepare fresh 1% agar LB plates with ampicillin. Take a colony from the “contaminated” TG1 control plate as well as the glycerol stock used to grow the log-phase TG1 and streak it on the new plates. If no growth is observed, this indicates that the ampicillin stock is old and that the TG1 was not infected with phage prior to this.

Problem 3: Only a single hit is found on the anti-M13 ELISA plate.

Possible cause: Only a few binders are present after the rounds of panning.

Solution: Another 93 colonies can be picked to prepare a second master plate to test in an anti-M13 ELISA.

Problem 4: None or very few anti-M13 ELISA hits present after the second round of panning using 5 nM of target protein.

Possible cause: Might not have enrichment.

Solution: Go back to 50 nM phage and pan using 10 nM antigen. If there are still not hits, try **Troubleshooting 3**. If there are still no hits, consider looking at another library. Use the seroconversion ELISA data to guide this

decision.

## Acknowledgments

This work was supported by the Rosalind Franklin Institute, with funding delivery partner the Engineering and Physical Sciences Research Council UK (EPSRC) and grants from the Biotechnology and Biological Sciences Research Council UK (BBSRC) for nanobody discovery (ref. BB/V018523/1) and the Wellcome Trust for technology development (ref. 223733/Z/21/Z). Graphical overview created with [BioRender](#).

## Competing interests

The authors declare that they have no competing interests with respect to the work described.

## Ethical considerations

Immunisations and handling of the llamas were performed under the authority of the UK Home Office project license PA1FB163A.

## References

1. Hamers-Casterman, C., Atarhouch, T., Muyldermans, S., Robinson, G., Hammers, C., Songa, E. B., Bendahman, N. and Hammers, R. (1993). [Naturally occurring antibodies devoid of light chains](#). *Nature* 363(6428): 446–448.
2. Bao, G., Tang, M., Zhao, J. and Zhu, X. (2021). [Nanobody: a promising toolkit for molecular imaging and disease therapy](#). *EJNMMI Res.* 11(1): 6.
3. Muyldermans, S. (2020). [A guide to: generation and design of nanobodies](#). *FEBS J.* 288(7): 2084–2102.
4. Parker, J. L., Deme, J. C., Wu, Z., Kuteyi, G., Huo, J., Owens, R. J., Biggin, P. C., Lea, S. M. and Newstead, S. (2021). [Cryo-EM structure of PepT2 reveals structural basis for proton-coupled peptide and prodrug transport in mammals](#). *Sci. Adv.* 7(35): eabh3355.
5. Girt, G. C., Lakshminarayanan, A., Huo, J., Dormon, J., Norman, C., Afrough, B., Harding, A., James, W., Owens, R. J., Naismith, J. H., et al. (2021). [The use of nanobodies in a sensitive ELISA test for SARS-CoV-2 Spike 1 protein](#). *R. Soc. Open Sci.* 8(9): e211016.
6. Huo, J., Mikolajek, H., Le Bas, A., Clark, J. J., Sharma, P., Kipar, A., Dormon, J., Norman, C., Weckener, M., Clare, D. K., et al. (2021). [A potent SARS-CoV-2 neutralising nanobody shows therapeutic efficacy in the Syrian golden hamster model of COVID-19](#). *Nat. Commun.* 12(1): 5469.
7. Akkermans, O., Delloye-Bourgeois, C., Peregrina, C., Carrasquero-Ordaz, M., Kokolaki, M., Berbeira-Santana, M., Chavent, M., Reynaud, F., Raj, R., Agirre, J., et al. (2022). [GPC3-Unc5 receptor complex structure and role in cell migration](#). *Cell* 185(21): 3931–3949.e26.
8. Pardon, E., Laeremans, T., Triest, S., Rasmussen, S. G. F., Wohlkönig, A., Ruf, A., Muyldermans, S., Hol, W. G. J., Kobilka, B. K., Steyaert, J., et al. (2014). [A general protocol for the generation of Nanobodies for structural biology](#). *Nat. Protoc.* 9(3): 674–693.
9. Kirchhofer, A., Helma, J., Schmidhals, K., Frauer, C., Cui, S., Karcher, A., Pellis, M., Muyldermans, S., Casas-Delucchi, C. S., Cardoso, M. C., et al. (2010). [Modulation of protein properties in living cells using nanobodies](#). *Nat. Struct. Mol. Biol.* 17(1): 133–138.
10. Maier, J., Traenkle, B. and Rothbauer, U. (2015). [Real-time analysis of epithelial-mesenchymal transition using fluorescent single-domain antibodies](#). *Sci. Rep.* 5(1): 13402.
11. Tonikian, R., Zhang, Y., Boone, C. and Sidhu, S. S. (2007). [Identifying specificity profiles for peptide](#)

- [recognition modules from phage-displayed peptide libraries](#). *Nat. Protoc.* 2(6): 1368–1386.
12. Le Bas, A., Mikolajek, H., Huo, J., Norman, C., Dormon, J., Naismith, J. and Owens, R. (2022). [Production and Crystallization of Nanobodies in Complex with the Receptor Binding Domain of the SARS-CoV-2 Spike Protein](#). *Bio Protoc.* 12(9): e4406.